

# DINAMAP<sup>®</sup> PRO Series 100-400V2 Monitor Operation Manual



**GE Medical Systems**  
Information Technologies

*We bring good things to life.*  
[gemedicalsystems.com](http://gemedicalsystems.com)



**DINAMAP® PRO 400V2 Monitor**

# DINAMAP® PRO Monitor

## Operation Manual

This manual is for DINAMAP® PRO Monitor Models 100V2, 200V2, 300V2, and 400V2, all with printers.

- PRO 100V2: BP and Pulse
- PRO 200V2: BP, Pulse, and Temp
- PRO 300V2: BP, Pulse, and SpO<sub>2</sub>
- PRO 400V2: BP, Pulse, Temp, and SpO<sub>2</sub>

The model of the Monitor determines which menu option buttons appear on the LCD. Please refer to applicable sections.

The Model Number 100-400V2 is generic in nature and reflects the range of product codes available. Your product may be labeled with a specific product code such as DINAMAP® PRO Model 410. "V2" refers to the second version of the product's design.

### ***Reissues and Updates***

Changes occurring between issues are addressed through Change Information Sheets, Addendums, and replacement pages. If a Change Information Sheet does not accompany this manual, it is correct as printed.

### ***Errors and Omissions***

If errors or omissions are found in this manual, please notify:  
GE Medical Systems *Information Technologies*  
4502 Woodland Corporate Boulevard  
Tampa, FL 33614  
1-877-274-8456

Part No. 2009802-001B

The content of this document including all figures and drawings is proprietary information of GE Medical Systems *Information Technologies*, provided solely for purposes of operation, maintenance or repair, and dissemination for

other purposes or copying thereof is prohibited without prior written consent by GE Medical Systems *Information Technologies*, Tampa, Florida.

Illustrations may show design models; production units may incorporate changes.

### ***Hierarchy of Warnings and Cautions***

A **general warning** is a statement that alerts the user to the possibility of injury, death, or other serious adverse reactions associated with the misuse of the device. A **warning** relates to steps in a procedure.

A **general caution** is a statement that alerts the user to the possibility of a problem with the device associated with its use or misuse. Such problems include device malfunction, device failure, damage to the device or damage to other property. A **caution** relates to steps in a procedure.

© GE Medical Systems *Information Technologies*, 2002. All rights reserved. Tampa, FL 33614

Printed in the U.S.A.

### **United States**

GE Medical Systems *Information Technologies*  
4502 Woodland Corporate Boulevard  
Tampa, FL 33614

### **European Representative**

GE Medical Systems  
Information Technologies GmbH  
Postfach 60 02 65  
D-79032 Freiburg, Germany  
Tel. +49 761 45 43 - 0  
Fax +49 761 45 43 - 233

# Contents

<b>Introduction .....</b>	<b>7</b>
About the DINAMAP® PRO Monitor .....	7
Product Compliance .....	10
Symbols .....	11
<b>Getting Started.....</b>	<b>13</b>
Unpacking the Monitor and Accessories .....	13
Power Sources .....	13
Powering the Monitor .....	13
Rear Panel Connections .....	16
Front Panel Controls and Indicators .....	17
Switching the Monitor On and Off.....	19
Liquid Crystal Display (LCD).....	19
Using the Printer .....	20
Installing the Paper .....	20
Printer Alarms .....	21
Storage .....	22
<b>Using the Monitor .....</b>	<b>23</b>
Noninvasive Blood Pressure Determination .....	23
Procedures.....	27
Manual Mode .....	29
Auto Mode.....	29
Stat Mode .....	30
TURBO*TEMP.....	31
Masimo SET® SpO <sub>2</sub> .....	37
NELLCOR® Oximax™ SpO <sub>2</sub> .....	46
<b>Using the Menu System .....</b>	<b>55</b>
Introduction.....	55
Liquid Crystal Display.....	55
Rotor .....	58
Menu Tree .....	58
Main Menu.....	58
Vitals Button (UK: All Obs) .....	59
More... Button .....	60
Set BP Button (UK: BP Mode).....	60
Alarms Button .....	62
Trend Button .....	63
Print Button .....	65
More... Menu .....	66
SpO <sub>2</sub> Button .....	66
Config Button.....	67
Pwr Sav (Sleep Mode) .....	68
Time .....	69
Rotor .....	69
Temp.....	70
Display Button .....	70
Service Button.....	71
Clinician Menu .....	72
Error and Warning Messages .....	75
Alarms Button .....	75
OK Button.....	76

<b>Appendix A</b> .....	<b>77</b>
Technical Specifications.....	77
BP .....	77
Temperature.....	78
SpO <sub>2</sub> .....	78
Mechanical.....	82
Power Requirements.....	82
Environmental .....	83
<b>Appendix B</b> .....	<b>85</b>
Patient Alarms.....	85
System Alarms .....	85
Failsafe Alarm.....	85
Hierarchy of Alarms.....	86
<b>Appendix C</b> .....	<b>91</b>
Principles of Noninvasive Blood Pressure Determination.....	91
<b>Appendix D</b> .....	<b>95</b>
Reorder Codes.....	95
<b>Appendix E</b> .....	<b>97</b>
Warranty, Service, and Spare Parts.....	97
Repairs.....	98
Packing Instructions.....	98
Service Manuals .....	98
<b>Appendix F</b> .....	<b>99</b>
Maintenance .....	99
Cleaning the Monitor .....	99
Cuff Cleaning and Disinfection.....	99
Temperature Devices.....	100
SpO <sub>2</sub> Sensors.....	100
Storage and Battery Care .....	101
Fuses .....	102
Calibration .....	102
Leak Testing.....	102
Disposal of Product Waste .....	103
<b>Appendix G</b> .....	<b>105</b>
Host Port Connector (rear panel).....	105

# Introduction

## ***About the DINAMAP® PRO Monitor***

DINAMAP® PRO Monitors provide noninvasive determination of systolic blood pressure, diastolic blood pressure, mean arterial pressure, pulse rate, temperature, and oxygen saturation. These portable AC- and DC-operated monitors are primarily intended for use in hospital acute care settings such as outpatient surgery, accident and emergency, labor and delivery, GI/endoscopy, and medical/surgical units.

The PRO Monitor comes in four different models: PRO 100V2, 200V2, 300V2, and 400V2, all with printers.

- PRO 100V2: BP and Pulse
- PRO 200V2: BP, Pulse, and Temp
- PRO 300V2: BP, Pulse, and SpO<sub>2</sub>
- PRO 400V2: BP, Pulse, Temp, and SpO<sub>2</sub>

All of the main operations of the PRO Monitor are easy to use. Please review the factory default settings and, where applicable, enter settings appropriate for your use. The “Using the Monitor” section of this manual explains how to use the Monitor in its most simple form, while the “Using the Menu System” section explains how to customize measurements by using the menu system.

## ***Indications***

The PRO Monitor is intended to monitor one patient at the bedside.

## ***Contraindications***

This device is not designed, sold, or intended for use except as indicated.

Federal law (U.S.A.) restricts this device to sale by or on the order of a clinician.

## ***Warnings***

- **Do not use the PRO Monitor in the presence of magnetic resonance imaging (MRI) devices. There have**

been reports of sensors causing patient burns when operating in an MRI environment.

- Do not use the Monitor in the presence of flammable anesthetics.
- To help prevent unintended current return paths with the use of high frequency (HF) surgical equipment, ensure that the HF surgical neutral electrode is properly connected.
- To avoid personal injury, do not perform any servicing unless qualified to do so.
- **WARNING:** These Monitors should not be used on patients who are connected to cardiopulmonary bypass machines.
- If powering the Monitor from an external power adapter or converter, use only power adapters and converters approved by GE Medical Systems *Information Technologies*.
- The Monitor does not include any user-replaceable fuses. Refer servicing to qualified service personnel.
- To reduce the risk of electric shock, do not remove the cover or the back. Refer servicing to a qualified service person.
- If the accuracy of any determination reading is questionable, first check the patient's vital signs by alternate means and then check the PRO Monitor for proper functioning.

### *Cautions*

- Do not use replacement batteries other than the type supplied with the Monitor. Replacement batteries are available from GE Medical Systems *Information Technologies*. See Appendix D.
- The PRO Monitor is designed to conform to Electromagnetic Compatibility (EMC) standard IEC 601-1-2, 1993 and will operate accurately in conjunction with other medical equipment which also meets this requirement. To avoid interference problems affecting the Monitor, do not use the Monitor in the

presence of equipment which does not conform to these specifications.

- **Place the PRO Monitor on a rigid, secure surface. Monitor must only be used with mounting hardware, poles, and stands recommended by GE Medical Systems *Information Technologies*. See Appendix D.**
- **The weight of the accessory basket contents should not exceed 6.6 lb (3 kg).**
- **Arrange the power cord, air hoses, and all cables carefully so they do not constitute a hazard.**
- **Verify calibration of BP parameter (temp and pulse oximeter do not require calibration). Ensure that the display is functioning properly before operating the PRO Monitor.**
- **Do not immerse the Monitor in water. If the Monitor is splashed with water or becomes wet, wipe it immediately with a dry cloth.**
- **Do not gas sterilize or autoclave.**
- **The PRO Monitor, when used with applied parts and accessories *approved by GE Medical Systems Information Technologies*, is protected against defibrillator damage.**

## Notes

- Waveforms may be distorted and readings inaccurate when electrosurgical cautery equipment is used while monitoring with the PRO Monitor.
- The electromagnetic compatibility profile of the PRO Monitor may change if accessories other than those specified for use with the PRO Monitor are used.
- Trend data are retained in the PRO Monitor when it is turned off, except when the default is overridden by selecting the Trend button under the Service menu.

### ***Product Compliance***

The DINAMAP® PRO Monitor is classified in the following categories for compliance with IEC 601-1:

- Class I, internally powered
- Transportable
- For continuous operation
- Not suitable for use in the presence of flammable anesthetics
- Not for use in the presence of an oxygen-enriched atmosphere (oxygen tent)
- Type BF applied parts
- IPX1, degree of protection against ingress of water
- Sterilization/Disinfection, see Appendix F



DINAMAP® PRO MONITOR  
CLASSIFIED WITH RESPECT TO ELECTRIC SHOCK, FIRE  
AND MECHANICAL AND OTHER SPECIFIED HAZARDS  
ONLY IN ACCORDANCE WITH CAN/CSA C22.2 NO.  
601.1. ALSO EVALUATED TO IEC-601-2-30.



**0086**


This product conforms with the essential requirements of the Medical Device Directive. Accessories without the CE mark are not guaranteed to meet the Essential Requirements of the Medical Device Directive.

## Symbols

The following symbols are associated with the PRO Monitor.

**Note:** The type of model determines which symbols appear on the Monitor.

	Attention, consult accompanying documents
	Defibrillator-proof type BF equipment
	Power ON/OFF
	SpO <sub>2</sub>
	START/STOP BP
	AUTO BP
	Battery Power
	MAP (Mean Arterial Pressure)
	Temperature
	Beats Per Minute
	
	SILENCE
	External Communications Port Connector
	External AC or DC Power Indicator
	External DC Power Input
	External AC Power Input
	SpO <sub>2</sub> Motion Artifact / Low Perfusion (US)
	SpO <sub>2</sub> Motion Artifact / Low Perfusion (Europe)

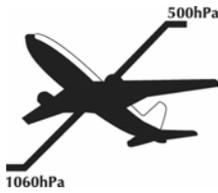

Packaging label depicting the transportation and storage atmospheric pressure range of 500 to 1060 hPa.

## IPX1

The DINAMAP® PRO Monitor is protected against vertically falling drops of water and conforms to IEC-529 standard at level of IPX1. Vertically falling drops of water shall have no harmful effects to the Monitor.

# Getting Started

## ***Unpacking the Monitor and Accessories***

Before attempting to use the PRO Monitor, take a few minutes to become acquainted with the Monitor and its accessories. Unpack the items carefully, and check them against the contents checklist enclosed in one of the accessory boxes. This is also a good time to check for any damage or shortage. If there is a problem or shortage, contact GE Medical Systems *Information Technologies*.

It is recommended that all the packaging be retained, in case the Monitor must be returned for service in the future.

## ***Power Sources***

The PRO Monitor is designed to operate from either an internal lead-acid battery, AC mains or an IEC 601-1 compliant DC power source (see Appendix A). For replacement rechargeable batteries, please refer to the Service section of this manual.

The Monitor contains five fuses. Two AC line input fuses are mounted internally and are replaceable only by qualified service personnel. The remaining three fuses are auto-resettable and mounted within the Monitor. These fuses protect the low voltage DC input, the battery, and the +5 V output on the host port connector.

## ***Powering the Monitor***

Before the PRO Monitor is used for the first time, the battery should be charged in the Monitor for at least 8 hours.

Refer to the illustration of the rear panel connections. Looking at the rear of the PRO Monitor, remove the battery compartment cover. Insert the rechargeable battery into the compartment so that the battery terminals fit into the power clips at the bottom of the compartment. Then replace the cover. Insert the power cord plug into the mains external power socket (2) and plug into an AC outlet.

Refer to the illustration of the front panel controls and indicators. With mains or external DC power connected, the green external power indicator LED (14) will light to indicate that external power is being applied and that the battery is

charging. If the battery is not inserted, the external power indicator LED will flash (short flash approx. every 4 sec). When the Monitor is running on battery power, a battery icon appears in LCD area 3 (toggling with the time indicator) indicating the charge status.

During battery-only operation, the yellow battery power indicator LED (17) will light. When the battery becomes discharged beyond the low battery threshold, the indicator will begin to flash, and the Monitor will sound warning beeps every 30 seconds. At this point, the Monitor should be connected to an AC outlet to recharge the battery. If the Monitor continues to be used without charging the battery, the message **WARNING: THE BATTERY IS TOO LOW FOR MONITOR TO FUNCTION. TURN MONITOR OFF** appears. The Monitor shuts down all functions until it is turned off and the battery is recharged or replaced. To run the Monitor on AC power, it must be powered off and then on again.

Battery charging will take place as long as the Monitor remains connected to an external AC power source. A battery that is fully discharged can be fully recharged in 1 hour 50 minutes when the Monitor is switched off or 8 hours if the Monitor is switched on.

### Notes

- To prolong the life of the battery, keep the Monitor connected to an AC outlet whenever possible. NEVER allow the battery to become completely discharged. A fully charged battery will power the Monitor for approximately 2 hours and should survive between 200 and 500 charge/discharge cycles. When it is necessary to replace the battery, refer to the “Compatibility Table and Reorder Codes” listed in Appendix D. To ensure full charge cycles, replace only with a recommended battery. If the Monitor is to be stored for some time, first charge the battery and then remove it and store it separately from the Monitor.
- For continued safety, use only a power cord of listed type SJT, three-conductor, min. No. 18 AWG, terminated


# Getting Started

in a medical/hospital grade attachment plug, provided with the following cord tag: "Hospital Grade Plug." Grounding integrity can only be maintained when equipment is connected to an equivalent receptacle marked "Hospital Grade."

- Where the integrity of the external earth conductor in the installation or its arrangement is in doubt, the Monitor must be operated from its internal battery.

## *General Caution*


- **Do not touch either the pin of the DC input connector (3) or the terminals within the battery compartment (1) and the patient at the same time.**

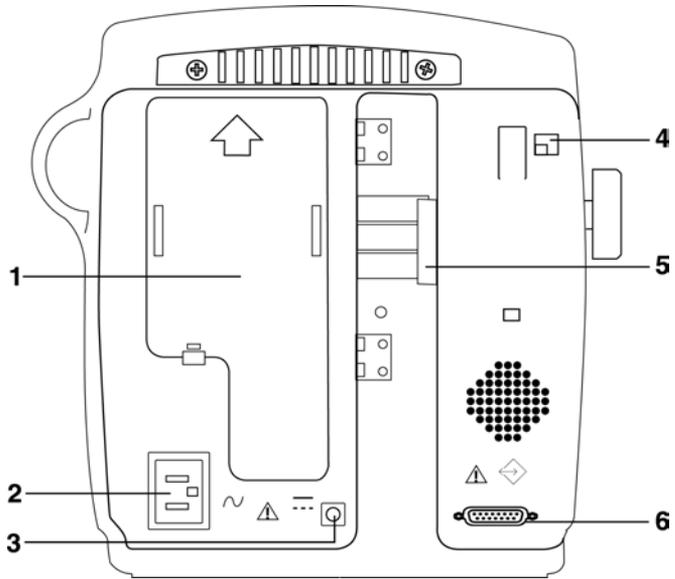

### ***Rear Panel Connections***

- 1** Battery compartment cover: Retains and protects internal battery.
- 2** Mains input: Used to connect to AC power supply.
- 3** External power socket: To be used with approved GE Medical Systems *Information Technologies* AC-DC power converter ONLY.
- 4** Inactive temperature cable storage: Inactive temperature probe cable attaches here (Models 200V2 and 400V2).
- 5** Pole clamp: Used to clamp monitor to pole or stand.
- 6** Data interface connector: Host communications port (15 way D-type RS-232 serial port) for use only with equipment conforming to IEC 601-1, configured to comply with IEC 601-1-1.

# Getting Started

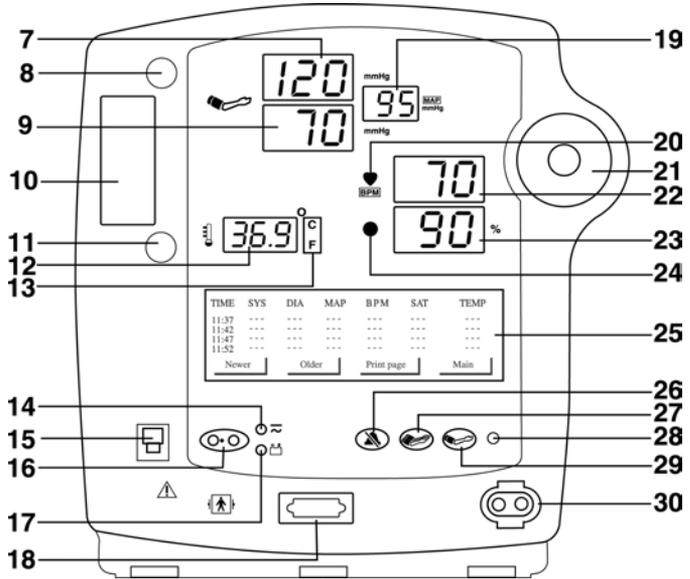

## Front Panel Controls and Indicators

- 7 Systolic pressure display: 3-digit red LED indicates measured systolic BP in mmHg.
- 8 Active temperature probe holster: Temperature probe that is being used stored here (Models 200V2 and 400V2).
- 9 Diastolic pressure display: 3-digit red LED indicates measured diastolic BP in mmHg.
- 10 Temperature probe cover storage: Box of probe covers stored here (Models 200V2 and 400V2).
- 11 Inactive temperature probe holster: Extra temperature probe can be stored here (Models 200V2 and 400V2).
- 12 Temperature display: 4-digit red LED indicates measured temperature (Models 200V2 and 400V2).
- 13 °C / °F display: Indicates whether temperature is being displayed in degrees Celsius or Fahrenheit (Models 200V2 and 400V2).

- 14** External power indicator: Green LED indicates external power status and battery charging status of monitor.
- 15** Temperature probe connector: Temperature probe cable attaches here (Models 200V2 and 400V2).
- 16** ON/OFF switch: Controls on/off state of monitor; push for power on and push again for power off.
- 17** Battery power indicator: Yellow LED indicates operation and charge status of internal battery.
- 18** SpO<sub>2</sub> sensor connector: SpO<sub>2</sub> sensor extension cable attaches here (Models 300V2 and 400V2).
- 19** Mean arterial pressure display: 3-digit red LED indicates measured MAP in mmHg and shows instantaneous cuff pressure during BP determination.
- 20** SpO<sub>2</sub> pulse indicator: Yellow LED in heart symbol flashes to indicate that real-time pulse rate measurements are being derived from SpO<sub>2</sub> signals (Models 300V2 and 400V2).
- 21** Rotor: Used to highlight and select items in LCD menus; if monitor is off, pressing rotor will switch monitor on.
- 22** Pulse BPM display: 3-digit yellow LED shows pulse rate in beats per minute.
- 23** SpO<sub>2</sub> display: 3-digit red LED indicates oxygen saturation in % (Models 300V2 and 400V2).
- 24** SpO<sub>2</sub> motion/artifact indicator LED: For NELLCOR, LED illuminates when motion artifact is detected (Models 300V2 and 400V2). For Masimo, LED illuminates when low perfusion or low signal quality is detected. (Models 300V2 and 400V2).
- 25** LCD (liquid crystal display): Displays all alarms, user interface messages, and configuration options.
- 26** Alarm silence switch: Alternately mutes and enables audible alarms; when pushed once after alarm sounds (silence on), switch lights to indicate that audible alarms have been silenced for 2 minutes.
- 27** AUTO BP key: Press to start Auto BP mode.
- 28** Light sensor: Automatically measures ambient light to set LED display intensity.
- 29** START/STOP BP key: Press to start or stop a BP, Auto, Stat, or Vitals determination.
- 30** Cuff connector: BP cuff hose attaches here.

# Getting Started

## *Switching the Monitor On and Off*




To switch the DINAMAP PRO Monitor on, press and hold the power ON/OFF switch (16) for at least 10 seconds or press the rotor (21).

As the Monitor powers up, it will run a short self-test routine, which will flash all the indicator lights and then beep the warning speaker. After a few seconds the system will be ready for operation, as indicated by the appearance of the main menu on the LCD (25).

To switch the Monitor off, push the power ON/OFF switch (16) again. This will terminate any measurements that may be in progress and automatically deflate the cuff.

When the Monitor is operating on the internal battery only, battery life is enhanced by the use of the sleep mode. However, the PRO Monitor will not enter sleep mode if an alarm is active. If no controls are used and no determinations are being made, the Monitor will enter sleep mode after a time which can be preset by the operator. All LED displays will be blanked except for a dash in the far-left systolic position, and any existing readings will be transferred to the LCD, which displays the message "Sleep Mode Active." Moving the rotor or pressing a key will "wake up" the Monitor.

## *Liquid Crystal Display (LCD)*

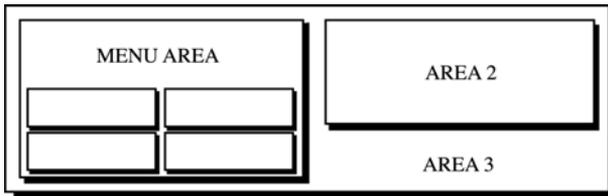

### **Menu Area**

This area displays the name of the menu that has option buttons available for selection. Normal text in the menu area

appears dark on a light background, while the text of selected buttons appears light on a dark background.

**Note:** Some menus have six option buttons. In these cases, there is no space available to display the menu title.

### **Area 2**

This area displays data from one of three different sources.

- Source 1: SpO<sub>2</sub> plethysmograph (Models 300V2 and 400V2)
- Source 2: Last three BP readings
- Source 3: Error and warning messages

**Note:** Refer to “Display Button” in the “Using the Menu System” section for instructions on setting Area 2.

### **Area 3**

This area displays the time, the time lapsed since the last BP determination, the battery icon (if operating on battery power, the time and battery icon toggle), and the BP and Printer modes.


## ***Using the Printer***

### **Installing the Paper**

Turn the PRO Monitor so that the side is facing you. While grasping the side of the Monitor, firmly press the notched indentations on the printer door to open it. The printer door will pop open. With the Monitor powered on, place the roll of paper into the compartment so that the end of the paper comes off the top, and thread it between the two printer plates. As the paper touches the plates, the paper will begin to auto-feed itself into the printer. Feeding the end of the paper strip through the slot in the door, firmly press the notched indentation on the side of the printer door to close it. Use the paper release lever to clear a paper jam or manually feed the paper.

**Note:** Make sure that the roll of paper is tightly wound.

# Getting Started

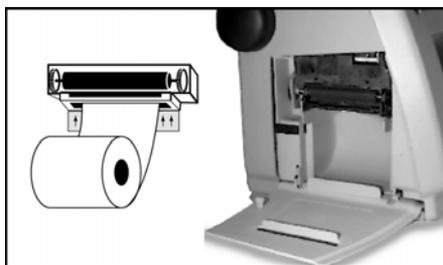

Any time the paper is loaded, the printer automatically prints a test strip with the DINAMAP® PRO name on it. If no print is visible on the paper, check that the paper roll has been installed in the correct position (refer to diagram). To tear off the printout, use a slight sideways action to pull the paper sharply up across the serrated edge of the door.

## Printer Alarms

If the Monitor is switched on with no paper installed or with the printer door open, the message “No Paper” will appear next to “PRNT” in Area 3 of the LCD. When new paper is installed and the printer door is closed, the message will change to “Manual” for Manual print or “Auto” for Auto print, depending on the status before the paper change.

If the paper runs out during a print request or if an attempt is made to print when no paper is installed, the message “Printer - No Paper” will appear in Area 2 of the LCD and an audible alarm will sound. In addition, the message “No Paper” will appear next to “PRNT” in Area 3 of the LCD. To clear the alarm, press the rotor. The message in Area 3 of the LCD will remain until new paper is installed and the printer door is closed. (See “Using the Menu System.”)

Installing new paper will cause the DINAMAP PRO header to be printed, thereby confirming that the paper is installed correctly and that the printer is operational. The message next to “PRNT” in Area 3 of the LCD will change to “Auto” or “Manual” to identify the operating mode of the printer. After power-off, the operating mode of the printer returns to the previous user-selected setting (Auto or Manual) unless specified otherwise in the Print button under the Service Button.

### **Storage**

Store thermal paper in a cool, dry place. The printed strip (thermal paper recording) should not be

- exposed to direct sunlight,
- exposed to temperatures over 100 °F/38 °C or relative humidity over 80%, or
- placed in contact with adhesives, adhesive tapes, or plasticizers such as those found in all PVC page protectors.

**Note:** When in doubt about long-term storage conditions, store a photocopy of the thermal paper recording.

### ***Cautions***

- **The paper is thermally activated; therefore, do not store it in a hot place as discoloration may result.**
- **Use only replacement paper rolls (58 mm) from GE Medical Systems *Information Technologies*.**

# Using the Monitor

## *Noninvasive Blood Pressure Determination*

### **Description**

The BP parameter is included in Models 100V2, 200V2, 300V2, and 400V2. Blood pressure is monitored noninvasively in the PRO Monitor by the oscillometric method, which measures the amplitude of the pressure oscillations within the blood pressure cuff. Further information about the oscillometric method is in Appendix C.

The PRO Monitor has four BP modes: 1. Manual, 2. Auto, 3. Stat, and 4. Vitals (UK: All Obs). The mode, which is selected by the user, is shown on the LCD (25). The BP measurements are automatic, and once the cycle is complete the LED displays (7, 9, 19, 22) show systolic pressure, diastolic pressure, mean arterial pressure, and pulse rate.

1. Manual BP determinations are started by pressing the START/STOP BP key (29). In the Manual mode, the blood pressure is determined one time.
2. Auto BP determinations are started by selecting the AUTO BP key (27) or the Auto button under the Set BP (UK: BP Mode) button in the Main menu.

When Auto mode is selected, a number at the right of the Auto button indicates the time interval between each reading. To change the time interval, choose the box around the number and turn the rotor until the desired interval is reached. The interval can be set between 1 and 120 minutes (1, 2, 3, 4, 5, 10, 15, 20, 30, 45, 60, 90, and 120 minutes). Press the rotor to confirm the setting.

3. Stat determinations are started by selecting the Stat button under the Set BP button (UK: BP Mode) in the Main menu. In the Stat mode, the blood pressure is determined as many times as possible in 5 minutes.
4. Vitals (UK: All Obs) determinations are started by selecting the Vitals (UK: All Obs) button in the Main menu. (Refer to the "Using the Menu System" section.) Selection of this button initiates a BP determination while

allowing SpO<sub>2</sub> and predictive temperature determinations to be monitored and recorded (depending on Monitor model). In the Vitals (UK: All Obs) mode, the blood pressure is determined one time.

Before each BP determination, the Monitor performs a test to ensure that the cuff pressure is below a specified level. The determination is delayed until this condition is met. During the delay, the BP values are displayed as zero.

The Monitor senses the type of hose being used and automatically uses adult/pediatric monitoring parameters or neonatal monitoring parameters, as appropriate.

Audible and visible alarms occur when a value for systolic pressure, diastolic pressure, mean arterial pressure, or pulse rate is outside the selected high or low limit.

Instructions for cleaning and disinfecting BP cuffs are in Appendix F.

### ***General Warnings***

- **The PRO Monitor will not measure blood pressure effectively on patients who are experiencing seizures or tremors.**
- **Arrhythmias will increase the time required by the PRO Monitor to determine a blood pressure and may extend the time beyond the capabilities of the Monitor.**
- **In Manual mode, the PRO Monitor displays the results of the last blood pressure determination for a duration of time set by the user or until another determination is completed. If a patient's condition changes between one determination and the next, the Monitor will not detect the change or indicate an alarm condition.**
- **Devices that exert pressure on tissue have been associated with purpura, skin avulsion, compartmental syndrome, ischemia and/or neuropathy. To minimize these potential problems, especially when monitoring at frequent intervals or over extended periods of time, make sure the cuff is applied appropriately and**

# Using the Monitor

examine the cuff site and the limb distal to the cuff regularly for signs of impeded blood flow.

- Do not apply external pressure against cuff while monitoring. Doing so may cause inaccurate blood pressure values.
- Use care when placing cuff on extremity used to monitor other patient parameters.
- The PRO Monitor is designed for use only with dual-tube cuffs.
- Use only accessories recommended by GE Medical Systems *Information Technologies*. Failure to use recommended accessories may result in inaccurate readings. See Appendix D.
- Blood pressure cuffs should be removed from the patient when the Monitor is powered off. If the extremity remains cuffed under these conditions or if the interval between blood pressure determinations is prolonged, the patient's limb should be observed frequently and the cuff placement site should be rotated as needed.

## *General Cautions*

- Accuracy of BP measurement depends on using a cuff of the proper size. It is essential to measure the circumference of the limb and to select the proper size cuff. The air hoses are color-coded according to size of the patient. The gray 12- or 24-foot hose (3.66 m or 7.3 m) is required on patients who require cuff sizes from infant through thigh cuffs. The teal (blue-green) 12-foot hose (3.66 m) is required for the neonatal cuff sizes #1 through #5.
- If it becomes necessary to move the cuff to another limb, make sure the appropriate size cuff is used.
- The pulse rate derived from a BP determination may differ from the heart rate derived from an EKG waveform because the PRO Monitor measures actual peripheral pulses, not electrical signals or contractions from the heart. Differences may occur because electrical signals at the heart occasionally fail to

**produce a peripheral pulse or the patient may have poor peripheral perfusion. Also, if a patient's beat-to-beat pulse amplitude varies significantly (e.g., because of pulsus alternans, atrial fibrillation, or the use of a rapid-cycling artificial ventilator), blood pressure and pulse rate readings can be erratic, and an alternate measuring method should be used for confirmation.**

### **General Notes**

- A patient's vital signs may vary dramatically during the use of cardiovascular agents such as those that raise or lower blood pressure or those that increase or decrease heart rate.
- Because treatment protocols based on the patient's blood pressure may rely on specific values and differing measurement methods, such as auscultatory, clinicians should note a possible variance from values obtained with the PRO Monitor in planning patient care management. The PRO Monitor values are based on the oscillometric method of noninvasive blood pressure measurement and correspond to comparisons with intra-aortic values within ANSI /AAMI Standards for accuracy (a mean difference of  $\pm 5$  mmHg, and a standard deviation of  $\pm 8$  mmHg).
- Several conditions may cause the BP parameter to calculate and display only the mean arterial pressure (MAP) without a systolic and diastolic reading. These conditions include very low systolic and amplitude fluctuations, so an accurate calculation for these values can't be made (e.g., patient in shock); too small of a difference between systolic and MAP calculations in relationship to the difference between diastolic and MAP; or a leak has occurred in the PRO Monitor (1. Check all BP connections 2. Monitor may need calibration and leak testing). If only the MAP value is displayed, the systolic and diastolic will display dashes (—) and an alarm message "N99-BP FAILED" will be displayed.

# Using the Monitor


## Procedures

1. Connect the end of the air hose which has quick-release clips to the cuff connector (30) on the front of the Monitor. Make sure that the hose is not kinked or compressed.


**Note:** To disconnect the hose from the Monitor, squeeze the quick-release clips together and pull the plug from the cuff connector (30).

2. Select the appropriate blood pressure measurement site. Because normative values are generally based on this site and as a matter of convenience, the upper arm is preferred. When upper arm size or shape, the patient's clinical condition, or other factors prohibit use of the upper arm, the clinician must plan patient care accordingly, taking into account the patient's cardiovascular status and the effect of an alternative site on blood pressure values, proper cuff size, and comfort. The figure shows the recommended sites for placing cuffs.

**Warning: Do not place the cuff on a limb being used for intravenous infusion or any area where circulation is compromised or has the potential to be compromised.**

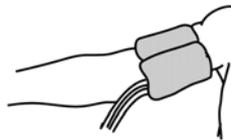

Adult/Pediatric

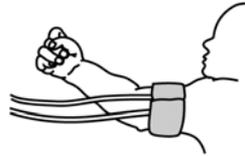

Neonate

3. If patient is standing, sitting, or inclined, ensure that cuffed limb is supported to maintain cuff at level of patient's heart. If cuff is not at heart level, the difference in systolic and diastolic values due to hydrostatic effect must be considered. Add 1.80 mmHg to values for every inch (2.54 cm) above heart level. Subtract 1.80 mmHg from values for every inch (2.54 cm) below heart level.
4. Select appropriate cuff size. Measure patient's limb and select appropriately sized cuff according to size marked on cuff or cuff packaging. When cuff sizes overlap for a specified circumference, choose the larger size cuff.

**Precaution: Accuracy depends on use of proper size cuff.**

5. **Inspect cuff for damage. Replace cuff when aging, tearing, or weak closure is apparent. Do not inflate cuff when unwrapped.**

**Precaution: Do not use cuff if structural integrity is suspect.**

6. Connect the cuff to the air hose. Thread the cuff connectors onto the hose connectors until finger tight. Do not overtighten.

**Warning: It is mandatory that the appropriate hose and cuff combination be used. Any attempt to modify the hose will inhibit the Monitor from switching between the neonatal and adult measurement modes.**

**Note:** In normal use, each cuff will have its own hose, so it will not usually be necessary to disconnect them. If it is necessary to do so, carefully unscrew the cuff from the hose. Care should be taken in reconnecting the cuff to a hose, ensuring that threads of the cuff and hose are in alignment and no cross-threading occurs.

7. Inspect patient's limb prior to application.

**Precaution:** Do not apply cuff to areas where skin is not intact or tissue is injured.

8. Palpate artery and place cuff so that patient's artery is aligned with cuff arrow marked "artery."
9. Squeeze all air from cuff and confirm that connection is secure and unoccluded and that tubing is not kinked.
10. Wrap cuff snugly around the patient's limb. Cuff index line must fall within the range markings. Ensure that hook and loop closures are properly engaged so that pressure is evenly distributed throughout cuff. If upper arm is used, place cuff as far proximally as possible.
11. Proper cuff wrapping should be snug, but should still allow space for a finger between patient and cuff. Cuff should not be so tight as to prevent venous return between determinations.  
**Warning: Using a cuff that is too tight will cause venous congestion and discoloration of the limb, but using a cuff that is too loose may result in no readings and/or inaccurate readings.**
12. Proceed with monitoring in the Manual, Auto, Stat, or Vitals (UK: All Obs) mode.

# Using the Monitor

## Manual Mode



To start a determination, press the START/STOP BP key (29). A normal, uninterrupted Manual cycle takes about 40 seconds. The cuff pressure must drop below 5 mmHg (neonate) or 15 mmHg (adult) before another determination can be started. BP information will be displayed for a duration of time set by the user. The factory default setting is 2 minutes (refer to “Retain Values” in the *Using the Menu System* section for information on how to change this setting). The user can set the display duration in 2, 5, 10, 15, 20, 30, 45, 60, 90, 120 minute increments. If another determination is initiated during the user-set timeframe, the display will change. This applies to Manual and Vitals (UK: All Obs) modes. After power-off, the operating mode returns to the default setting of Manual. The default setting of Manual can be overridden to return to the previous user-selected setting (Auto or Manual) by selecting Set BP (UK: BP Mode) under the Service menu.

**Note:** The START/STOP BP key is an on-off switch; pressing it will stop any BP determination (Manual, Auto, Stat, or Vitals) that is in progress.

## Auto Mode



Auto BP determinations are started by selecting the AUTO BP key (27) or the Auto button under the Set BP button (UK: BP Mode) in the Main menu.

Selecting the AUTO BP key (27) brings up the Set BP menu (UK: BP Mode) and automatically starts an Auto BP determination as long as the Monitor is in Manual BP mode. If the Monitor is already in Auto BP mode, selection of the AUTO BP key (27) brings up the Set BP menu (UK: BP Mode) without starting a new determination until the preset time interval has expired. Pressing the START/STOP BP key during a series of Auto BP determinations will cancel the determination in progress.

When Auto mode is selected, a number at the right of the Auto button indicates the time interval between each reading. To change the time interval, choose the box around the number and turn the rotor until the desired interval is reached. The interval can be set between 1 and 120 minutes (1, 2, 3, 4, 5, 10, 15, 20, 30, 45, 60, 90, and 120 minutes). Press the rotor to confirm the setting. After power-off, the operating mode returns to the default setting of Manual. The default setting of Manual can be overridden to return to the previous user-selected setting (Auto or Manual) by selecting Set BP (UK: BP Mode) under the Service menu.

In the Auto mode, the pressure must be below 5 mmHg (neonate) or 15 mmHg (adult) for at least 30 seconds before another determination can be started. BP information will be displayed on the LED until the next determination is started. This applies to Auto mode only.

**Note:** To cancel an Auto BP determination, select the Manual button in the Set BP menu (UK: BP Mode).

### **Stat Mode**

Multiple BP readings can be taken at any time by selecting the Stat button under the Set BP button (UK: BP Mode) in the Main menu. Stat mode can also be accessed by pressing the AUTO BP key (27) and then selecting the Stat button when the Set BP menu (UK: BP Mode) appears.

If a Manual determination is not in progress, a 5-minute series of determinations will start. If a Manual determination is in progress, that determination will become the first in the series. A normal, uninterrupted Stat sequence will give the first set of systolic, diastolic, and mean arterial pressure values and pulse rate within 15 to 20 seconds. Selecting the Stat button during a series of Stat determinations will cancel the determination in progress and the rest of the series. BP information will be displayed on the LED until the determination has been canceled or completed. This applies to Stat mode only.

The series begins with cuff inflation to a pressure above the previous systolic pressure or, if no previous systolic value is stored, to approximately 150 mmHg for adult/pediatrics. The

# Using the Monitor

initial target pressure selection for neonates is 110 mmHg. Artifact rejection is relaxed in the Stat mode for adult/pediatric patients to allow for accelerated determinations. If a BP or Stat reading has been made previously, the first new systolic value will flash on the LED display (7) within a few seconds and will continue to flash until the end of the determination. At that point a short tone will sound and the updated systolic, diastolic, and mean arterial pressures and pulse rate will appear on their LED displays (7, 9, 19, 22). The Monitor will begin another determination once the pressure is below 5 mmHg for 8 seconds (neonates) or 15 mmHg for 4 seconds (adults), unless the 5-minute period has ended or the determination has been canceled.

**Note:** Alarm limits are disabled while in Stat mode.

## **TURBO★TEMP**



### **Description**

The Temperature parameter is included in Models 200V2 and 400V2. The PRO Monitor uses IVAC\* TURBO★TEMP™ technology and can be used with both oral and rectal temperature probes. The TURBO★TEMP™ parameter consists of an electronic thermometer that uses a heat-sensing device known as a thermistor to sense temperature. The thermistor is part of the electrical circuit and is located at the tip of the probe. To obtain temperatures, the probe tip measures the rate of change in temperature when the thermistor comes into contact with surrounding tissue. A final temperature is calculated based on this rate of change.

During a temperature determination, the temperature display (12) provides a progress meter and probe ready indicator. In the far-left position, a single horizontal line indicates the probe is ready to start a determination after removal from the probe holster. In the far-right position of the temperature display, a “chase sequence” around the outside space indicates a predictive temperature determination is in progress. During monitor mode, the temperature readings flash constantly.

\*IVAC is a trademark of Alaris Medical Systems

Temperature is shown on the temperature display in degrees Celsius or Fahrenheit, and the unit of measure is indicated by the °C °F display (13). The default, which is Celsius, can be changed in the Clinician Menu (please refer to the “Using the Menu System” section of this manual).

Temperature determinations can be either predictive or monitoring. The user can make this selection from the Main Menu using the More... and then Temp buttons (please refer to the “Using the Menu System” section of this manual).

### **Predictive Mode**

In predictive mode, a final temperature is displayed with an audible single tone. Upon initiation of a measurement, the old temperature measurement, if present, is cleared. A predictive measurement is initiated when the probe is removed from the probe holster. When a predictive mode measurement is in progress, it is terminated when one of the following occurs:

- A final value is determined.
- The probe is inserted into the probe holster.
- The temperature measurement mode is automatically switched to monitor mode because a predictive result was unable to be computed.
- A Temperature alarm is issued.

A predictive temperature measurement value is automatically cleared when it is older than the display time set by the user and the probe has been returned to the probe holster. When the value is older than 5 minutes and the probe has not been returned to the probe holster, that value is clear.

### **Monitor Mode**

Monitor mode is most commonly used for axillary temperature determinations. In monitor mode, the display is updated continually as the patient’s temperature rises or falls. Monitor mode is automatically initiated when a predictive mode measurement terminates after approximately 10 seconds of not being able to successfully compute a result.

# Using the Monitor

When a monitor mode measurement is in progress, it is terminated when the probe is inserted into the probe holster. Note: These temperature readings are not stored in trends and not reported via host comms.

## ***General Warning***

- **The performance of the Monitor may be degraded if it is operated outside of the environmental conditions specified in Appendix A.**


## ***General Cautions***

- **Be careful not to overextend the coiled cord of the temperature probe. Overextension can damage the probe coil connector interfaces.**
- **Accurate oral temperatures (blue) can only be obtained by placing the probe under the tongue in the right or left sublingual pocket. Temperatures in other locations in the mouth can vary by more than 2 °F or 1 °C.**
- **Accurate rectal temperatures can only be obtained by using the red temperature probe. Red and blue temperature probes are *not* interchangeable.**
- **Do not allow the tip of the predictive temperature probe to come into contact with a heat source (e.g., hands or fingers) prior to taking a temperature determination. If this occurs, allow 5 seconds for the probe tip to cool before proceeding.**
- **Electromagnetic Compatibility: Operating the thermometer near equipment which radiates high-energy electromagnetic and radio frequencies (electrosurgical/cauterizing equipment, portable radios, cellular telephones, etc.) may cause false alarm conditions. If this happens, reposition the thermometer away from the source of interference and perform a new measurement.**
- **Use only IVAC® Probes and P850A probe covers with the TURBO★TEMP™ Thermometer. The size, shape, and thermal characteristics of the probe covers can affect the performance of the instrument. Inaccurate**



# Using the Monitor

7. Remove the probe. Discard the disposable cover by holding the probe as you would a syringe and pressing the button on the probe handle. Place the probe in the probe holster. An audible single-tone sounds. Once you place the probe in the probe holster, the temperature values will be cleared in the user selected timeframe.

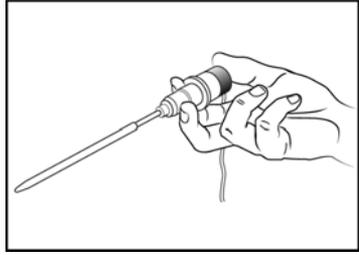

## Procedures for Rectal Predictive Mode Determinations

To measure rectal temperature use the red rectal probe.

1. Connect the temperature probe cable to the temperature probe connector.
2. Remove the temperature probe from the probe holster. An audible single-tone sounds. Place a protective temperature probe cover on the probe. Failure to firmly install the probe cover may result in the probe cover becoming loose or disengaged during use. Be careful not to press the probe ejection button where the cord exits the probe as this might loosen or eject the probe cover.
3. Touch the tissue about a half inch (1.3 cm) above the sphincter muscle and carefully insert the probe, using current hospital technique for penetration. (The use of a lubricant is optional.)
4. The determination begins automatically. To ensure continuous tissue contact and maximize patient comfort, hold the probe in position until the determination is complete. This takes approximately 10 seconds, during which time a "chase sequence" of arrows in the right-side of the **Temperature** window appears to indicate progress.
5. When the determination is complete, an audible tone sounds and the temperature appears on the display.
6. Remove the probe. Discard the disposable cover by holding the probe as you would a syringe and pressing the button on the probe handle. Place the probe in the probe holster. An audible single-tone sounds. Once you place the probe in the probe holster, the temperature values will be cleared in the user selected timeframe.

## Procedures for Monitor Mode Determinations (Axillary Determinations)

1. Connect the temperature probe cable to the temperature probe connector.
2. Remove the temperature probe from the probe holster. An audible single-tone sounds. Place a protective temperature probe cover on the probe and insert the probe in the patient's axilla, making sure the tip of the probe is in contact with the skin and positioned as close as possible to the axillary artery with the patient's arm held close to his/her side.
5. The determination begins automatically. The dashes in the vital signs area are replaced with a flashing temperature value as soon as the probe warms up. Leave the probe in place for the same length of time as required by standard hospital procedure for taking an axillary temperature. The Monitor does not beep to indicate a final reading.
6. Remove the probe. Discard the disposable cover by holding the probe as you would a syringe and pressing the button on the probe handle. Place the probe in the probe holster. An audible single-tone sounds. Once you place the probe in the probe holster, the temperature values will be cleared in the user-selected timeframe.



### Notes

- If there is a long delay from the time the probe is removed from the probe holster until it is inserted into the patient's mouth, it is possible that the instrument will not display a final temperature. If this occurs, insert the probe into the probe holster, remove it again, and start a new measurement.
- If an alarm is actively sounding, an audible tone will not sound.
- If the probe is not inserted/placed within 30 seconds of probe removal, the Monitor may issue an alarm.

# Using the Monitor

- If the probe tip temperature is 92.0° F or higher (33.3° C) when taken out of the probe holster, the thermometer will not be able to perform a predictive measurement. Instead, the thermometer will automatically go into monitor mode. The temperature reading will then flash. A correct final temperature reading may require 3 minutes or longer. The Monitor will not beep at final temperature. It will continue to monitor the patient's temperature until the probe is removed from the patient and returned to the probe holster.

## *Masimo Set\* SpO<sub>2</sub>*



### **Description**

The SpO<sub>2</sub> parameter in the PRO Monitor is available in two different leading technologies: NELLCOR® OxiMax™ and MASIMO SET®. Please refer to the front of your PRO Monitor to see which SpO<sub>2</sub> technology your monitor contains. The SpO<sub>2</sub> technology logo will be on the front fascia. **This section refers to MASIMO SET® SpO<sub>2</sub> technology.**


The SpO<sub>2</sub> parameter is included in Models 300V2 and 400V2. To begin SpO<sub>2</sub> monitoring, simply place the SpO<sub>2</sub> sensor on the patient's finger; monitoring begins automatically. Functional oxygen saturation (SpO<sub>2</sub>) of arterial blood is noninvasively and continuously monitored in the PRO Monitor using pulse oximetry technology from MASIMO SET®. Functional SpO<sub>2</sub> is the ratio of oxygenated hemoglobin to hemoglobin that is capable of transporting oxygen. This ratio, expressed as a percentage, is shown in the SpO<sub>2</sub> window, which is continually updated.

\*Masimo SET® is a trademark of Masimo Corporation. Possession or purchase of this device does not convey any express or implied license to use the device with replacement parts which would, alone, or in combination with this device, fall within the scope of one or more of the patents relating to the device.

Pulse rate derived from SpO<sub>2</sub> appears in the Pulse Rate window, and the SpO<sub>2</sub> pulse indicator flashes synchronization with the real-time pulse rate measurements that are derived from the SpO<sub>2</sub> signal. A tone sounds at a rate corresponding to the pulse rate and at a pitch corresponding to the SpO<sub>2</sub> saturation level. The pitch is highest at 100% oxygen saturation, and it becomes lower as the saturation level falls. The Monitor displays a pulse bar graph. The pulse bar graph is proportional to the signal quality. The artifact indicator LED (24) illuminates when low perfusion or poor signal quality is detected.

Audible and visible alarms occur when SpO<sub>2</sub> levels are outside the alarm limits. When a limit alarm occurs, an alarm message appears in SpO<sub>2</sub> window.


**Note:** Limit alarms, printing, and trending are not available for the first 10 seconds of SpO<sub>2</sub> monitoring.

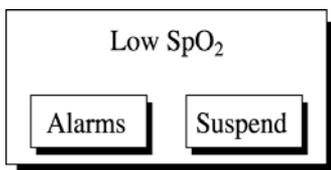

If you select the Alarms button, the Alarms menu appears. This menu is used to adjust the violation limits for BP and SpO<sub>2</sub>. Refer to "Alarms Button" in the "Using the Menu System" section.

If you select the Suspend button, the SpO<sub>2</sub> alarm is suspended for 2 minutes and then the PRO Monitor returns to normal SpO<sub>2</sub> monitoring. A message informing the user that SpO<sub>2</sub> is suspending appears in Area 2 and dashes appear in the SpO<sub>2</sub> LED while the SpO<sub>2</sub> alarm suspend is counting down. Selecting **Cancel** will cancel the SpO<sub>2</sub> alarm suspension and return to monitoring SpO<sub>2</sub>.

# Using the Monitor

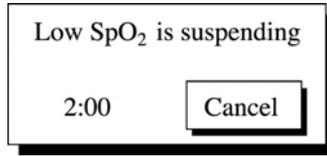

If the Monitor is unable to detect a pulse for 10 seconds during normal SpO<sub>2</sub> monitoring, the values in the LED flash, alternating patient values with dashes. The Monitor returns to normal SpO<sub>2</sub> reporting of values when several consecutive good pulse determinations are made.

## *General Cautions*

- **If any measurement does not seem reasonable, first check the patient's vital signs by alternate means and the check the pulse oximeter for proper functioning.**
- **Inaccurate measurements may be caused by:**
  - **Incorrect sensor application or use**
  - **Significant levels of dysfunctional hemoglobins (e.g. carboxyhemoglobin or methemoglobin)**
  - **Intravascular dyes such as indocyanine green or methylene blue.**
  - **Exposure to excessive illumination, such as surgical lamps (especially ones with xenon light source), bilirubin lamps, fluorescent lights, infrared heating lamps, or direct sunlight (exposure to excessive illumination can be corrected by covering the sensor with a dark or opaque material)**
  - **excessive patient movement**
  - **venous pulsations**
  - **placement of sensor on an extremity with a blood pressure cuff, arterial catheter, or intravascular line**
- **The pulse oximeter can be used during defibrillation, but the readings may be inaccurate for a short time.**
- **loss of pulse signal can occur in any of the following situation:**
  - **the sensor is too tight**

- there is excessive illumination from light sources such as surgical lamp, a bilirubin lamp, or sunlight
- a blood pressure cuff is inflated on the same extremity as the one with a SpO<sub>2</sub> sensor attached.
- the patient has hypotension, severe vasoconstriction, severe anemia, or hypothermia
- there is arterial occlusion proximal to the sensor
- the patient is in cardiac arrest or is in shock

### Sensors

- Before use, carefully read the Masimo® sensor directions for use.
- Use only Masimo® oximetry sensors for SpO<sub>2</sub> measurements. Other oxygen transducers (sensors) may cause improper SpO<sub>2</sub> performance.
- Tissue damage can be caused by incorrect application or use of an Masimo® sensor, for example by wrapping the sensor too tightly. Inspect the sensor site as directed in the sensor Directions for Use to ensure skin integrity and correct positioning and adhesion of the sensor.
- Do not use damaged Masimo® sensors. Do not use an Masimo® sensor with exposed optical components. Do not immerse the sensor in water, solvents, or cleaning solutions (the sensors and connectors are not waterproof). Do not sterilize by irradiation, steam, or ethylene oxide. See the cleaning instructions in the directions for use for reusable Masimo® sensors.
- Do not use damaged patient cables. Do not immerse the patient cables in water, solvents, or cleaning solutions (the patient cable connectors are not waterproof). Do not sterilize by irradiation, steam, or ethylene oxide. See the cleaning instructions in the directions for use for reusable Masimo® patient cables.
- Do not use the SpO<sub>2</sub> function during magnetic resonance imaging (MRI). Adverse reactions include potential burns to patients as a result of contact with attachments heated by the MRI radio frequency pulse, potential degradation of the magnetic resonance image, and potential reduced accuracy of SpO<sub>2</sub>

# Using the Monitor

measurements. Always remove oximetry devices and attachments from the MRI environment before scanning a patient.

- The use of cardio-green and other intravascular dyes at certain concentrations may affect the accuracy of the SpO<sub>2</sub> measurement.
- The SpO<sub>2</sub> function is calibrated to read functional arterial oxygen saturation. Significant levels of dysfunctional hemoglobins such as carboxyhemoglobin or methemoglobin may affect the accuracy of the SpO<sub>2</sub> measurement.
- The use of accessories, transducers, and cables other than those specified may result in increased emission and/or decreased immunity and inaccurate readings of the Monitor.
- As with any clip-on sensor, pressure is exerted. The clinician should be cautious in using a clip-on sensor on patients with compromised circulation (e.g., because of peripheral vascular disease or vasoconstricting medications).
- Do not perform any testing or maintenance on a sensor while it is being used to monitor a patient.
- Bright light sources (e.g., infrared heat lamps, bilirubin lights, direct sunlight, operating room lights) may interfere with the performance of the SpO<sub>2</sub> function. To prevent such interference, cover the sensor with opaque material.

## General Notes

- A patient's vital signs may vary dramatically during the use of cardiovascular agents such as those that raise or lower blood pressure or those that increase or decrease pulse rate.
- The PRO Monitor that is labeled with MASIMO SET<sup>®</sup> Technology is compatible only with MASIMO SET<sup>®</sup> sensors.
- Software development, software validation, and Risk and Hazard Analysis has been performed to a registered quality system.

- Use or patient-applied parts are latex-free.

### ***Procedures***

1. Select a sensor that is appropriate for the patient and the clinical situation.

**Warning: Do not use a damaged sensor or one with exposed electrical contacts.**

**Note:** Use only MASIMO<sup>®</sup> sensors, which are available from Masimo<sup>®</sup> Corporation and GE Medical Systems - *Accessories and Supplies*.

2. Following the directions for use supplied with the sensor, apply the sensor to the patient.

### ***Warnings***

#### **Patient safety:**

- **If you fail to apply the sensor properly, the patient's skin could be injured or the ability of the PRO Monitor to measure oxygen saturation could be compromised. For example, a clip-on sensor should never be taped shut. Taping the sensor could damage the patient's skin or impair the venous return, thus causing venous pulsation and inaccurate measurement of oxygen saturation.**
- **Excessive pressure from the sensor may cause necrosis of the skin.**

#### **Monitor performance:**

- **When an SpO<sub>2</sub> sensor is on a limb that has a blood pressure cuff, the SpO<sub>2</sub> data will not be valid when the cuff is inflated. If SpO<sub>2</sub> readings are required during the entire blood pressure determination, attach the SpO<sub>2</sub> sensor to the limb opposite the one with the blood pressure cuff.**
- **Remove nail polish and artificial nails. Placing a sensor on a polished or an artificial nail may affect accuracy.**

### ***Cautions***

#### **Patient safety:**

- **Do not place any clip-on sensor in a patient's mouth or on a patient's nose or toe.**
- **Do not place a clip-on finger sensor on a patient's thumb or across a child's foot or hand.**

# Using the Monitor

- **Observe the sensor site to assure adequate distal circulation. Sensor sites should be checked at least every 2 hours and rotated at least every 4 hours.**

## **Monitor performance:**

- **Placing a sensor distal to an arterial line may interfere with adequate arterial pulsation and compromise the measurement of SpO<sub>2</sub>.**
  - **Place the sensor so that the LEDs and the photodiode are opposite each other.**
3. Plug the SpO<sub>2</sub> sensor into the SpO<sub>2</sub> sensor extension cable. Then plug the SpO<sub>2</sub> sensor extension cable into the SpO<sub>2</sub> sensor connector.
  4. Proceed with monitoring. SpO<sub>2</sub> determinations run continuously and can run simultaneously with other measurements.

## **Troubleshooting Masimo SET® SpO<sub>2</sub> Parameter**

This section discusses potential difficulties and suggestions for resolving them. If the difficulty persists, contact a qualified service person or your local representative.

The service manual, which is for use by qualified service personnel provides additional troubleshooting information.

**PROBLEM:** The pulse amplitude bar indicates a pulse, but no oxygen saturation or pulse rate values appear on the screen.

### **CAUSE:**

- Excessive patient motion may be making it impossible for the SpO<sub>2</sub> function to find a pulse pattern.
- The sensor may be damaged.
- The patient's perfusion may be too low to allow the SpO<sub>2</sub> function to measure saturation and pulse rate.

### **SOLUTION:**

*Check the patient.*

- If possible, keep the patient still; check whether the SpO<sub>2</sub> sensor is applied securely and properly, and replace it if necessary; move the sensor to a new site; or

use a disposable adhesive sensor that may tolerate more motion.

- Replace the sensor.

**PROBLEM:** The SpO<sub>2</sub> value or the pulse rate changes rapidly; the pulse amplitude bar is erratic.

**CAUSE:**

- Excessive patient motion may be making it impossible for the SpO<sub>2</sub> function to find a pulse pattern.
- An electrosurgical unit (ESU) may be interfering with performance.

**SOLUTION:**

*Check the patient.*

- If possible, keep the patient still; check whether the sensor is applied securely and properly, and replace it if necessary; move the sensor to a new site; use a sensor that tolerates more motion.

If an ESU is interfering:

- Move the SpO<sub>2</sub> cable as far from the ESU as possible.
- Plug the Monitor and the ESU into different AC circuits.
- Move the ESU ground pad as close to the surgical site as possible.
- The sensor may be damp or may need to be replaced with a new sensor.

**PROBLEM:** The oxygen saturation measurement does not correlate with the value calculated from a blood gas determination.

**CAUSE:**

- The SpO<sub>2</sub> calculation may not have correctly adjusted for the effects of pH; temperature; CO<sub>2</sub>; fetal hemoglobin; or 2,3-DPG.
- Accuracy can be affected by incorrect sensor application or use; intravascular dyes; bright light; excessive patient movement; venous pulsations; electrosurgical interference; and placement of a sensor on an extremity that has a blood pressure cuff, arterial catheter, or intravascular line.

# Using the Monitor

## **SOLUTION:**

- Check that calculations have been corrected appropriately for the relevant variable. In general, calculated saturation values are not as reliable as direct laboratory hemoximeter measurements.
- If there is excessive light, cover the sensor with opaque material.
- Circulation distal to the sensor site should be checked routinely. The site must be inspected every 8 hours to ensure adhesion, skin integrity, and correct optical alignment. If skin integrity changes, move the sensor to another site.
- Try to keep the patient still, or change the sensor site to one with less motion.
- Observe all instructions, warnings, and cautions in this manual and in the directions for use of the sensor.

**PROBLEM:** A valid SpO<sub>2</sub> signal was present but has disappeared.

## **CAUSE:**

- A BP determination on the same limb is in progress.

## **SOLUTION:**

- An alarm message codes will appear on the screen, and the audible alarm will sound immediately.

**PROBLEM:** An SpO<sub>2</sub> PLACEMENT error has occurred.

## **CAUSE:**

- Weak or “noisy” signal.

## **SOLUTION:**

Check the patient.

- If possible, keep the patient still; check whether the sensor is applied securely and properly, and replace it if necessary; move the sensor to a new site; or replace the sensor.
- Change sensor type.
- Consider increasing perfusion using heat.

- If there is excessive light, cover the sensor with opaque material.

**PROBLEM:** A sensor error indicating a bad sensor has occurred.

**CAUSE:**

- The sensor or cable may be defective, or the cabling may be improperly connected.

**SOLUTION:**

Check the patient.

- If possible, keep the patient still; check whether the sensor/cable is applied securely and properly, and replace it if necessary.
- Disconnect and reconnect the sensor.

**NELLCOR® OxiMAX™ SpO<sub>2</sub>**



### **Description**

The SpO<sub>2</sub> parameter in the PRO Monitor is available in two different leading technologies: NELLCOR® OxiMAX and MASIMO SET®. Please refer to the front of your Monitor to see which SpO<sub>2</sub> technology your monitor contains. The SpO<sub>2</sub> technology logo will be on the front fascia of the Monitor. **This section refers to NELLCOR SpO<sub>2</sub> technology.**

Functional oxygen saturation (SpO<sub>2</sub>) of arterial blood is noninvasively and continuously monitored in the PRO Monitor using pulse oximetry technology from NELLCOR®. Functional SpO<sub>2</sub> is the ratio of oxygenated hemoglobin to hemoglobin that is capable of transporting oxygen.


Pulse rate derived from SpO<sub>2</sub> appears in the Pulse BPM display (22), and the SpO<sub>2</sub> pulse indicator (20) flashes synchronization with the real-time pulse rate measurements

\*NELLCOR® is a trademark of Nellcor Puritan Bennett

# Using the Monitor


that are derived from the SpO<sub>2</sub> signal. A tone sounds at a rate corresponding to the SpO<sub>2</sub> saturation level. The pitch is highest at 100% oxygen saturation and becomes a little lower with each 1% decrease in saturation level. The Monitor can display a pulse amplitude bar and a plethysmographic waveform on the LCD (25). The pulse amplitude bar graph is proportional to the arterial blood flow. The artifact indicator LED (24) lights continuously when the Monitor detects motion sufficient enough to affect readings.

Audible and visible alarms occur when SpO<sub>2</sub> levels are outside the alarm limits. When a limit alarm occurs, a message appears in Area 2 of the LCD display.

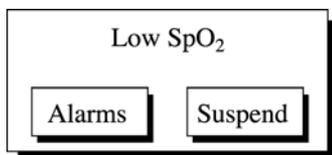

If you select the Alarms button, the Alarms menu appears. This menu is used to adjust the violation limits for BP and SpO<sub>2</sub>. Refer to "Alarms Button" in the "Using the Menu System" section.

If you select the Suspend button, the SpO<sub>2</sub> alarm is suspended for 2 minutes and then the PRO Monitor returns to normal SpO<sub>2</sub> monitoring. A message informing the user that SpO<sub>2</sub> is suspending appears in Area 2 and dashes appear in the SpO<sub>2</sub> LED while the SpO<sub>2</sub> alarm suspend is counting down. Selecting **Cancel** will cancel the SpO<sub>2</sub> alarm suspension and return to monitoring SpO<sub>2</sub>.

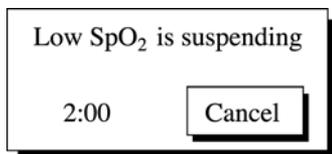

If the Monitor is unable to detect a pulse for 10 seconds during normal SpO<sub>2</sub> monitoring, the values in the LED flash, alternating patient values with dashes. The Monitor returns to

normal SpO<sub>2</sub> reporting of values when several consecutive good pulse determinations are made.

### ***General Warnings***

- **Do not use the SpO<sub>2</sub> function during magnetic resonance imaging (MRI). Adverse reactions include potential burns to patients as a result of contact with attachments heated by the MRI radio frequency pulse, potential degradation of the magnetic resonance image, and potential reduced accuracy of SpO<sub>2</sub> measurements. Always remove oximetry devices and attachments from the MRI environment before scanning a patient.**
- **The use of cardio-green and other intravascular dyes at certain concentrations may affect the accuracy of the SpO<sub>2</sub> measurement.**
- **The SpO<sub>2</sub> function is calibrated to read functional arterial oxygen saturation. Significant levels of dysfunctional hemoglobins such as carboxyhemoglobin or methemoglobin may affect the accuracy of the SpO<sub>2</sub> measurement.**
- **Pulse oximetry readings and pulse signals can be affected by certain environmental conditions, sensor application errors, and certain patient conditions. See the appropriate sections of this manual for specific safety information.**
- **The use of accessories, transducers, and cables other than those specified may result in increased emission and/or decreased immunity and inaccurate readings of the Monitor.**

### ***General Cautions***

- **As with any clip-on sensor, pressure is exerted. The clinician should be cautious in using a clip-on sensor on patients with compromised circulation (e.g., because of peripheral vascular disease or vasoconstricting medications).**

# Using the Monitor

- **Do not perform any testing or maintenance on a sensor while it is being used to monitor a patient.**
- **Bright light sources (e.g., infrared heat lamps, bilirubin lights, direct sunlight, operating room lights) may interfere with the performance of the SpO<sub>2</sub> function. To prevent such interference, cover the sensor with opaque material.**

## *General Notes*

- A patient's vital signs may vary dramatically during the use of cardiovascular agents such as those that raise or lower blood pressure or those that increase or decrease heart rate.
- When the front of the PRO Monitor is labeled with the Nellcor logo, use only NELLCOR sensors for SPO<sub>2</sub> measurements.
- Software development, software validation, and Risk and Hazard Analysis has been performed to a registered quality system.

## **Procedures**

1. Select a sensor that is appropriate for the patient and the clinical situation.

**Warning: Do not use a damaged sensor or one with exposed electrical contacts. Do not use sensors, cables, or connectors that appear damaged.**

**Note:** Use only NELLCOR<sup>®</sup> OxiMAX<sup>™</sup> sensors, which are available from: GE Medical Systems - *Accessories and Supplies*.

2. Following the directions for use supplied with the sensor, apply the sensor to the patient.

## **Warnings**

### **Patient safety:**

- **If you fail to apply the sensor properly, the patient's skin could be injured or the ability of the PRO Monitor to measure oxygen saturation could be compromised. For example, a clip-on sensor should never be taped shut. Taping the sensor could damage the patient's skin or impair the venous return, thus causing venous**

pulsation and inaccurate measurement of oxygen saturation.

- Excessive pressure from the sensor may cause necrosis of the skin.
- For additional warnings and information, refer to the NELLCOR<sup>®</sup> sensor's directions for use.

**Monitor performance:**

- When an SpO<sub>2</sub> sensor is on a limb that has a blood pressure cuff, the SpO<sub>2</sub> data will not be valid when the cuff is inflated. If SpO<sub>2</sub> readings are required during the entire blood pressure determination, attach the SpO<sub>2</sub> sensor to the limb opposite the one with the blood pressure cuff.
- Remove nail polish and artificial nails. Placing a sensor on a polished or an artificial nail may affect accuracy.

***Cautions***

**Patient safety:**

- Do not place any clip-on sensor in a patient's mouth or on a patient's nose or toe.
- Do not place a clip-on finger sensor on a patient's thumb or across a child's foot or hand.
- The sensor disconnect error message and associated alarm indicate that the sensor is either disconnected or the wiring is faulty. The user should check the sensor connection and, if necessary, replace the sensor, cable, or both.
- Observe the sensor site to assure adequate distal circulation. Sensor sites should be checked at least every 2 hours and rotated at least every 4 hours.

**Monitor performance:**

- Placing a sensor distal to an arterial line may interfere with adequate arterial pulsation and compromise the measurement of SpO<sub>2</sub>.
- Place the sensor so that the LEDs and the photodiode are opposite each other.

# Using the Monitor

3. Plug the SpO<sub>2</sub> sensor into the SpO<sub>2</sub> sensor extension cable. Then plug the SpO<sub>2</sub> sensor extension cable into the SpO<sub>2</sub> sensor connector (18).
4. Proceed with monitoring. SpO<sub>2</sub> determinations run continuously and can run simultaneously with other measurements.

## Troubleshooting the NELLCOR SpO<sub>2</sub> Parameter

This section discusses potential difficulties and suggestions for resolving them. If the difficulty persists, contact a qualified service person or your local representative.

The service manual, which is for use by qualified service personnel provides additional troubleshooting information.

**PROBLEM:** The pulse amplitude bar indicates a pulse, but no oxygen saturation or pulse rate values appear on the screen.

### CAUSE:

- Excessive patient motion may be making it impossible for the SpO<sub>2</sub> function to find a pulse pattern.
- The sensor may be damaged.
- The patient's perfusion may be too low to allow the SpO<sub>2</sub> function to measure saturation and pulse rate.

### SOLUTION:

*Check the patient.*

- If possible, keep the patient still; check whether the SpO<sub>2</sub> sensor is applied securely and properly, and replace it if necessary; move the sensor to a new site; or use a disposable adhesive sensor that may tolerate more motion.
- Replace the sensor.

**PROBLEM:** The SpO<sub>2</sub> value or the pulse rate changes rapidly; the pulse amplitude bar is erratic.

### CAUSE:

- Excessive patient motion may be making it impossible for the SpO<sub>2</sub> function to find a pulse pattern.
- An electrosurgical unit (ESU) may be interfering with performance.

**SOLUTION:**

*Check the patient.*

- If possible, keep the patient still; check whether the sensor is applied securely and properly, and replace it if necessary; move the sensor to a new site; use a sensor that tolerates more motion.

*If an ESU is interfering:*

- Move the SpO<sub>2</sub> cable as far from the ESU as possible.
- Plug the Monitor and the ESU into different AC circuits.
- Move the ESU ground pad as close to the surgical site as possible.
- The sensor may be damp or may need to be replaced with a new sensor.
- If the patient weighs less than 3 kg or more than 40 kg, apply an OXIMAX reusable sensor (except DS-100, OXI-A/N, OXI-P/I) or OxiCliq oxygen transducer to an appropriate site. These sensors have Faraday shields which provide added protection from high electronic noise and ambient light.

**PROBLEM:** The oxygen saturation measurement does not correlate with the value calculated from a blood gas determination.

**CAUSE:**

- The SpO<sub>2</sub> calculation may not have correctly adjusted for the effects of pH; temperature; CO<sub>2</sub>; fetal hemoglobin; or 2,3-DPG.
- Accuracy can be affected by incorrect sensor application or use; intravascular dyes; bright light; excessive patient movement; venous pulsations; electrosurgical interference; and placement of a sensor on an extremity that has a blood pressure cuff, arterial catheter, or intravascular line.

**SOLUTION:**

- Check that calculations have been corrected appropriately for the relevant variable. In general, calculated saturation values are not as reliable as direct laboratory hemoximeter measurements.

# Using the Monitor

- If there is excessive light, cover the sensor with opaque material.
- Circulation distal to the sensor site should be checked routinely. The site must be inspected every 8 hours to ensure adhesion, skin integrity, and correct optical alignment. If skin integrity changes, move the sensor to another site.
- Try to keep the patient still, or change the sensor site to one with less motion.
- Observe all instructions, warnings, and cautions in this manual and in the directions for use of the sensor.

**PROBLEM:** A valid SpO<sub>2</sub> signal was present but has disappeared.

**CAUSE:**

- A BP determination on the same limb is in progress.

**SOLUTION:**

- An alarm message (No signal) will appear on the screen, and the audible alarm will sound immediately.

**PROBLEM:** A bad signal error has occurred.

**CAUSE:**

- Weak or “noisy” signal.

**SOLUTION:**

*Check the patient.*

- If possible, keep the patient still; check whether the sensor is applied securely and properly, and replace it if necessary; move the sensor to a new site; or replace the sensor.
- Change sensor type.
- Consider increasing perfusion using heat.
- If there is excessive light, cover the sensor with opaque material.

**PROBLEM:** A sensor error indicating a bad sensor has occurred.

**CAUSE:**

- The sensor or cable may be defective, or the cabling may be improperly connected.

**SOLUTION:**

*Check the patient.*

- If possible, keep the patient still; check whether the sensor/cable is applied securely and properly, and replace it if necessary.
- Disconnect and reconnect the sensor.


# Using the Menu System

## ***Introduction***

The PRO Monitor is equipped with a liquid crystal display (25) and a rotor (21). Used together, these allow the operator to view and edit most of the Monitor's parameters and functions. When the Monitor is in use, a number of option buttons appear on the liquid crystal display (LCD). The model of the Monitor determines which menu option buttons appear on the LCD. The number of buttons and the specific options depend on the menu level. The rotor provides the means of choosing menu options and changing monitor settings.

## ***Liquid Crystal Display***

The LCD is divided into three areas, each of which has a distinct function.

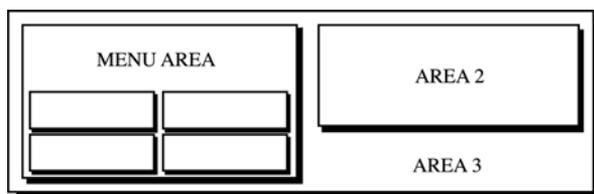

### **Menu Area**

This area displays the menu buttons that are available for selection. Normal text in the menu area appears dark on a light background, while the text of selected buttons appears light on a dark background.

### **Area 2**

This area displays BP and SpO<sub>2</sub> data and error and warning messages. The Display mode menu is used to select the data to be displayed.

### **Area 3**


This area displays the time, the time lapsed since the last Auto BP determination (if in Auto BP mode), the battery icon (if operating on battery power, the time and battery icon toggle), and the BP and printer modes.

**Note:** In cold ambient temperatures (below 50 °F / 10 °C), updates on the LCD can be delayed by approximately 1 second. This delay on the LCD does not affect the performance of the Monitor.



# Using the Menu System

## Clinician Mode Menu

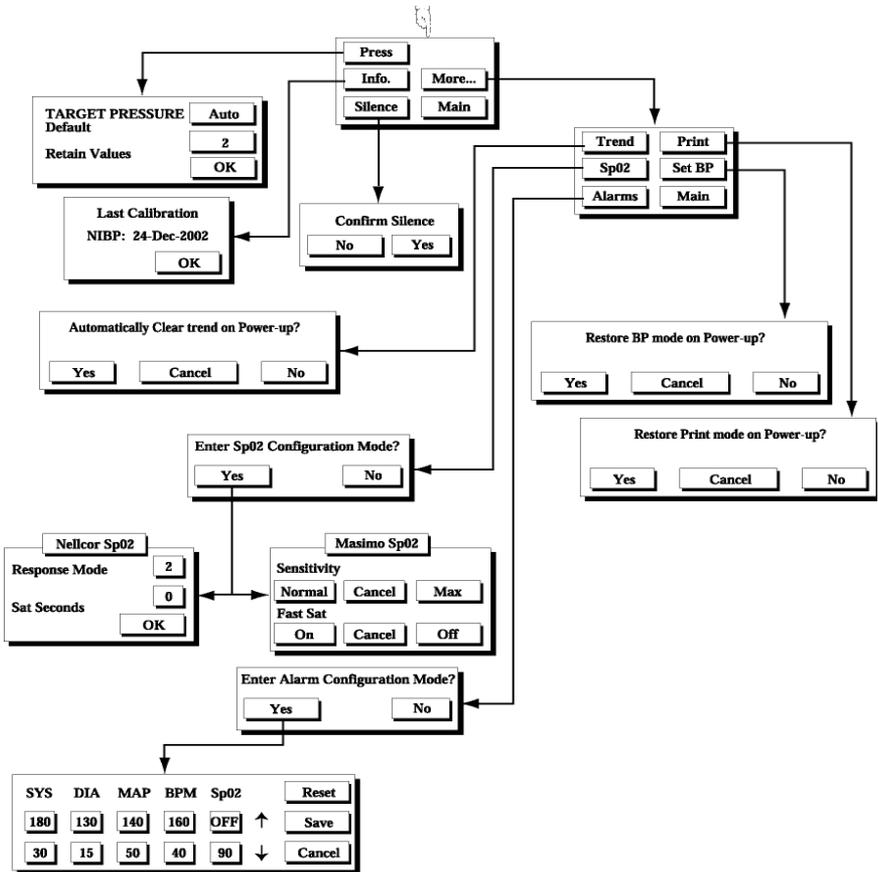

### ***Rotor***

Rotating the rotor causes option buttons to be highlighted (light text on a dark background). Turning the rotor produces a click. Turning it clockwise moves the highlighting clockwise over the available buttons, while turning it counterclockwise reverses the direction of the highlighting. Pressing the rotor selects the highlighted button and produces an audible tone.


Some menus (e.g., Alarms) contain values that can be changed by the operator. After the value is highlighted, the user selects it by pressing the rotor. Turning the rotor clockwise will cause the value to increase, and turning the rotor counterclockwise will cause the value to decrease. Pressing the rotor again will confirm the changed value.

### ***Menu Tree***

The menu tree on the previous page shows ***all*** possible choices available within the menu structure, from the top level downward.

### ***Main Menu***


This menu is the top level menu. It is displayed when the Monitor is first switched on and after the rotor has been inactive for 2 minutes, unless the Monitor is in sleep mode (Pwr Sav).

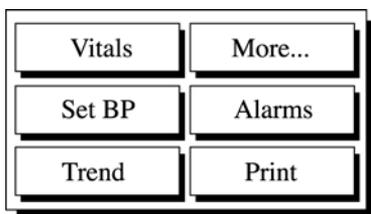

# Using the Menu System

## Vitals Button (UK: All Obs)

Selection of this button initiates a BP determination while allowing SpO<sub>2</sub> and predictive temperature determinations to be monitored and recorded (depending on Monitor model). When the Vitals determination is complete, a single “warble” sounds and all patient data are displayed on the LEDs and held for 2 minutes. The LCD shows\*:

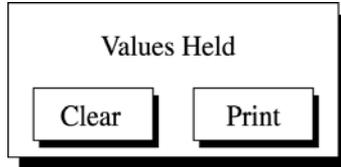

**Note:** If the printer is in “Auto” mode, the Print button does not appear as an option.

After this 2 minute static display, BP and temp values are displayed for the period of time pre-selected by the user for manual values (refer to “Retain Values” in the *Using the Menu System* section for information on how to change this setting). After the 2 minute static display, the SpO<sub>2</sub> parameter returns to active monitoring and dynamically displays the SpO<sub>2</sub> values.

## Notes

- If the Monitor is performing a Vitals determination, the Vitals button cannot be selected.
- If a BP determination is in progress, the Vitals button cannot be selected.
- A Vitals determination is canceled if the BP determination is canceled.
- A Vitals determination can be canceled by pressing either the AUTO BP or START/STOP key.
- During the 2 minute freeze period, SpO<sub>2</sub> monitoring and alarms are suspended.

### *Clear*

Selection of this button halts measurements and returns the user to the Main menu.

\*The model of the Monitor determines which menu option buttons appear on the LCD.

**Note:** If the SpO<sub>2</sub> plethysmograph is displayed on the LCD, the waveform pauses for 2 minutes or until the Clear button is selected. SpO<sub>2</sub> values are also retained the same manner as the BP and Temperature values.

### *Print*

Selection of this button causes the current data to be printed.

### **Notes**


- The Print button appears only when Print is set to Manual mode.
- If the printer is in Auto print mode, the data will be printed automatically.

### **More... Button**

Selection of this button displays the More... menu. The More... menu has six options (depending on model of Monitor), most of which have submenus. For this reason, instructions for the More... button are in a separate section.

### **Set BP Button (UK: BP Mode)**

Selection of this button displays the Auto, Stat, and Manual BP menu.

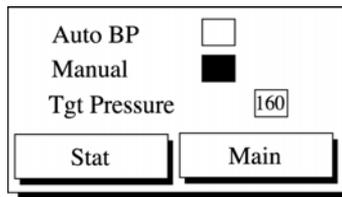

### *Auto*

Selection of this option starts an Auto BP determination. When Auto Mode is selected, a number at the right of the Auto button indicates the time interval between each reading. To change the time interval, choose the box around the number and turn the rotor until the desired interval is reached. The interval can be set between 1 and 120 minutes (1, 2, 3, 4, 5, 10, 15, 20, 30, 45, 60, 90, and 120 minutes). Press the rotor to confirm the setting. After power-off, the

# Using the Menu System

operating mode returns to the default setting of Manual. The default setting of Manual can be overridden to return to the previous user-selected setting (Auto or Manual) by selecting Set BP under the Service menu. To cancel an Auto BP determination, select the START/STOP key. To cancel an Auto BP mode, select the Manual button in the Set BP (UK: BP Mode) menu.

## *Manual*

Selection of this option starts a Manual BP determination. After power-off, the operating mode returns to the default setting of Manual. The default setting of Manual can be overridden to return to the previous user-selected setting (Auto or Manual) by selecting Set BP under the Service menu.

## *Tgt Pressure*

Selection of this option allows the user to set the BP target inflation pressure. The initial target pressure can be set between 100 and 250 mmHg in 5 mmHg increments. The factory default is 150 mmHg for adults and 110 for neonates. (This is indicated by “AUTO” at the end of adjustable range.) When using a neonate blood pressure cuff, if the target pressure is set to greater than 140 mmHg under the Set BP or Clinical menu, the Monitor automatically defaults to a target pressure of 110 mmHg. If the target pressure is set between 100 and 140 mmHg, then that setting is the target pressure that will be used. When the target pressure is changed, the next determination will use the new target inflation value if no systolic is available. Initial target pressure is restored to the factory default setting after power-off. The target pressure can be adjusted permanently in the Clinician menu of the Service mode (refer to “Press” in the “Using the Menu System” section).

## *Stat*

Selection of this option allows the user to start Stat determinations. When Stat is selected, blood pressure is determined as many times as possible in 5 minutes.

**Note:** Alarm limits are disabled while in Stat mode.

## *Main*

Selection of this button returns the user to the Main menu.

## Alarms Button

Selection of this button displays the Alarms menu. This menu is used to adjust the violation limits for BP, Pulse Rate, and SpO<sub>2</sub>. The values and ranges for these parameters are not stored when the Monitor is turned off. The user may edit the limits, but they are restored to the default values each time the Monitor is switched on. To permanently change the alarm limits, refer to “Alarms” under “Service Button” in the “Using the Menu System” section.


**Alarm Limits Table**

Parameter	Range	Default
Systolic High	35 - 245	180
Systolic Low	30 - 240	30
Diastolic High	15 - 195	130
Diastolic Low	10 - 190	15
MAP High	20 - 215	140
MAP Low	15 - 210	50
BPM High (Masimo)	35 - 235	160
BPM Low (Masimo)	30 - 230	40
BPM High (Nellcor)	35 - 250	150
BPM Low (Nellcor)	30 - 245	50
SpO <sub>2</sub> High	21 - 100	Off
SpO <sub>2</sub> Low	20 - 99	90

## Volume

Selection of this button displays the alarm volume submenu. The volume range is from 1 to 10, with 10 being the loudest. The alarm volume is stored when the Monitor is turned off and restored to the user’s preference each time the Monitor is switched on. Selection of the **Check** button allows the current volume setting to be heard. Selection of the **Main** button returns the user to the Main menu.

# Using the Menu System

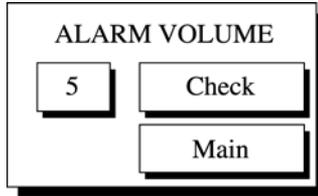

## *Auto*

Selection of this button updates the alarm limits relative to the current parameter values. Pressing this button cancels any limit violation alarm that is invalidated by this automatic limit change. The alarm limits are updated as follows:

Parameter Label	Label	High Limit	Low Limit
Systolic	SYS	SYS +30	SYS -30
Diastolic	DIA	DIA +30	DIA -30
MAP	MAP	MAP +30	MAP -30
Heart Rate	BPM	BPM +30	BPM -30
SpO <sub>2</sub>	SpO <sub>2</sub>	SpO <sub>2</sub> +5*	SpO <sub>2</sub> -5

\* If the reading plus the limit is greater than the valid range of measurement (e.g., SpO<sub>2</sub> +5 is greater than 100% ), the valid range of measurement becomes the limits.

## **Notes**


- In no case will the updated alarm limits be set beyond the valid limits in the Alarm Limits table.
- If no values are available, the limits will remain unchanged.

## *Main*

Selection of this button returns the user to the Main menu.

## **Trend Button**

Selection of this button displays the Trend mode menu.

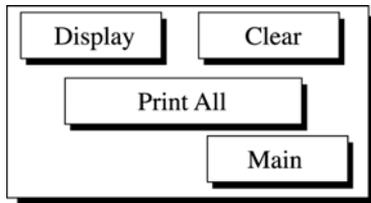

### *Display*

Selection of this button allows the operator to view the trend data.

**Note:** If the trend data have been lost (e.g., if the clock settings have been changed), the message “Trend Empty” will appear instead of the Newer, Older, and Print page buttons.

Newer and Older. These buttons may be used to move forward and backward through the recorded data. If no information is available, these buttons will not appear.

Print page. Selection of this button causes the displayed information to be printed. If no information is available, this button will not appear.


Main. Selection of this button returns the user to the Main menu.

### *Clear*

Selection of this button produces an advisory that the trend will be lost. Choosing **Yes** will erase the trend memory. Choosing **No** will retain the trend memory. This button disappears from the menu while printing and when Trend is empty.

### *Print All*

Selection of this button prints all the historical data available. When selected, this button temporarily changes to **Cancel** until the history has completed printing. Once printing is complete, the **Cancel** button returns back to the **Print All** button. This button disappears from the menu when Trend is empty.

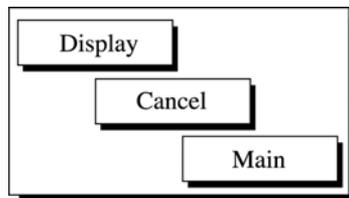


# Using the Menu System

## *Main*

Selection of this button returns the user to the Main menu.

## **Print Button**

Selection of this button displays the Print menu.

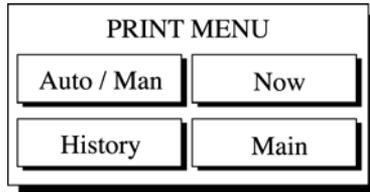

## *Auto/Man*

Pressing this button toggles between Automatic and Manual Printing modes. The current mode is displayed on Area 3 of the LCD. The Automatic mode prints the readings after each determination. The Manual mode, which is the factory default mode, requires the user to press the **Now** button to print the readings.

## *Now*

Selection of this button causes the last readings of the available parameters to be printed. If no readings are available, the message "No reading" is printed for that parameter. An error message appears if there is no paper in the printer.

## *History*


Selection of this button causes the entire contents of the trend memory to be printed. When selected, this button temporarily changes to **Cancel** until the history has completed printing.

## *Main*

Selection of this button returns the user to the Main menu.


### **More... Menu**

This menu is used to set the various operating modes of the Monitor.

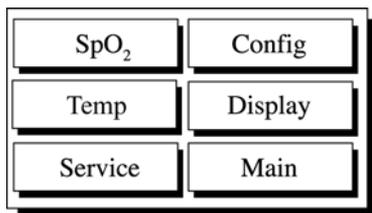

### **SpO<sub>2</sub> Button (Models 300V2 and 400V2)**

Selection of this button displays the SpO<sub>2</sub> mode menu, which is used to set the SpO<sub>2</sub> pulse tone volume and SpO<sub>2</sub> Average (Masimo ONLY).

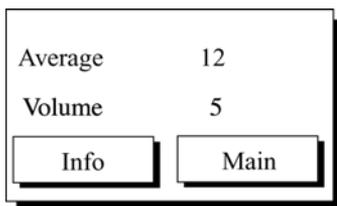

#### *Average*

For PRO Monitors with Masimo SET technology, the Average value is used in the calculation of SpO<sub>2</sub> values. Average can be set to 4, 6, 8, 10, 12, 14, or 16 seconds.

#### *Volume*

The pulse tone volume can be set in the range of **Off** to **9**. The value **Off** should be selected if no pulse tone is desired. The volume setting is stored when the Monitor is turned off and is restored to the user's preference each time the Monitor is switched on.

#### *Info*

With Nellcor SpO<sub>2</sub> technology installed, the info button displays the current settings for Response Mode and *SatSeconds™*.

# Using the Menu System

## *Response Mode (Nellcor Only)*

Response mode allows the user to specify the averaging technique to optimize measurements in the presence of various patient movement. Choose **Mode 1** (normal Response) when patients are active as in exercise protocols. Choose **Mode 2** (Fast response; default setting) for the general patient population.

## *SatSeconds™ (Nellcor Only)*

The *SatSeconds™* technique limit-controls the time that the % SpO<sub>2</sub> level may fall outside the alarm before an audible alarm sounds. Choose either **0, 10, 25, 50, or 100** seconds. If 0 is chosen this feature is disabled.

The *SatSeconds™* "Safety Net" is for patients with saturation levels having frequent excursions below the limit, but not staying below the limit long enough for the *SatSeconds* time setting to be reached. When 3 or more limit violations occur within 60 seconds, an alarm sounds even if the the *SatSeconds* time setting has not been reached.


With Masimo SET technology installed, the Info button displays the current settings for the Sensitivity and FastSAT features. Sensitivity can be set to Normal (the factory default) or Max. FastSAT can be set to On or Off (the factory default).

## *Main*

Selection of this button returns the user to the Main menu.


## **Config Button**

Selection of this button displays the Config mode menu, which allows the Power Save mode and time to be adjusted.

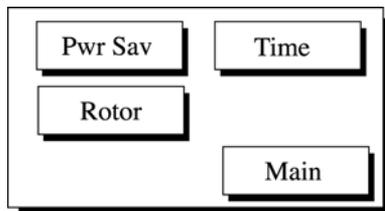

### *Pwr Sav (Sleep Mode)*

Selection of this button allows the operator to specify the time, in minutes, that elapses before the Monitor goes into “sleep” mode (LEDs blanked and LCD displaying values from LEDs). Sleep mode is available only if the Monitor is operating from its battery. Sleep mode conserves power while the Monitor is not in use. Once the Monitor is in Sleep mode, the user can return it to normal operation by touching any button or the rotor.

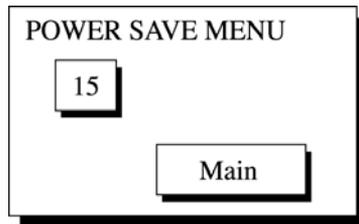

The Monitor enters Sleep mode only if the following are true:

- No alarm is active
- SpO<sub>2</sub> is not actively reporting patient statistics
- The keys and rotor have not been used for the preset time
- The Monitor is running from its battery
- No determinations are in progress
- The Monitor has been running from the battery for the entire preset time

The Monitor awakens from Sleep mode if any of the following occur:

- The rotor is turned or pressed
- Any of the keys are pressed
- An alarm condition is issued
- The battery supply level becomes discharged to a critical level
- A mains or suitable DC supply is connected
- An Auto BP or Temp determination starts
- An SpO<sub>2</sub> signal is detected
- A BP determination is started through the host comm

# Using the Menu System

## *Time*

Selection of this button allows the operator to change the internal time and date of the Monitor. The clock, which is maintained by an internal battery after power down, uses 24-hour format. The date is in the British format of dd/mmm/yyyy; however, to avoid confusion the month number has been substituted with a three-letter abbreviation. Leap years are calculated automatically.

The screenshot shows a menu interface for setting the time and date. At the top, there are two buttons: "Main" on the left and "Accept" on the right. In the center, the current time and date are displayed as "20 : 21 . 11 24-Aug-1999". Below this, there are six input fields arranged in two rows. The top row has labels "Hour", "Min", "Sec", "Day", "Mon", and "Year" above each field. The bottom row contains six empty rectangular boxes for entering the values.

Accept. Selection of this button produces an advisory to the user that the trend will be lost when the clock settings are changed. Choosing **Yes** will cause the Monitor to accept the new clock settings and erase the trend memory. Choosing **No** will cause the Monitor to retain the existing clock settings and the trend memory. Either choice returns the user to the Main menu.

Main. Selection of this button returns the user to the Main menu.

## *Rotor*

Selection of this button displays a panel for setting the volume of the beep that sounds when the rotor is turned. The range of adjustment is **Off** (default) to **9**, and the setting is retained when the Monitor is turned off.


The screenshot shows a menu interface for setting the rotor beep volume. The text "ROTOR BEEP" is centered at the top. Below it, the word "Volume" is on the left and the number "3" is on the right. At the bottom center, there is a button labeled "Main".

### *Main*

Selection of this button returns the user to the Main menu.

### **Temp**

Selection of this button displays the temperature submenu, which allows the user to choose the temperature label. When **C** (Celsius) is selected, the °C indicator lights. When **F** (Fahrenheit) is selected, the °F indicator lights.

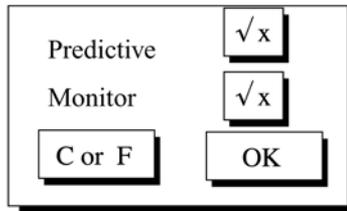

*C or F.* Selection of this button toggles the temperature display between Celsius and Fahrenheit.


*Predictive.* Selection of this button toggles between the two modes of temperature measurement. TURBO★TEMP predictive measurement initiates when the probe is removed from the holster and terminates when either a value is determined, the probe is re-holstered, or after a 2-minute time span when the software defaults to “Monitor” mode.

*Monitor.* Monitor mode is most commonly used for axillary temperature determinations. In monitor mode, the display is updated continually as the patient’s temperature rises or falls.

### **Display Button**

Selection of this button displays the Display mode menu. This menu is used to specify whether Area 2 of the LCD will display SpO<sub>2</sub> or BP data. If neither SpO<sub>2</sub> nor 3 NIBP is selected, Area 2 of the LCD will remain blank except for the pulse amplitude bar (if SpO<sub>2</sub> data are available) and any error or warning messages that may appear. The Display mode setting is maintained when the Monitor is switched off and on.

# Using the Menu System

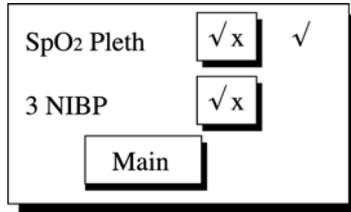

## *SpO<sub>2</sub> Pleth*

When this option is checked and SpO<sub>2</sub> data are available, the plethysmograph waveform and the pulse amplitude bar will be displayed.

## *3 NIBP*

When this option is checked, the last 3 NIBP readings will be displayed. If SpO<sub>2</sub> data is available, the pulse amplitude bar will also be displayed.

## *Main*


Selection of this button returns the user to the Main menu.

## **Service Button**


Selection of this button displays a keypad that allows the clinician to access some parts of the Service mode menu. To access the Clinician menu, use the rotor to select the numbers 1, 2, 3, 4 sequentially.

### **Notes**

- SpO<sub>2</sub> is automatically disabled when entering Service mode.
- Service modes that affect the calibration or alignment of the instrument are not available to the user. These modes are described in the Service Manual.

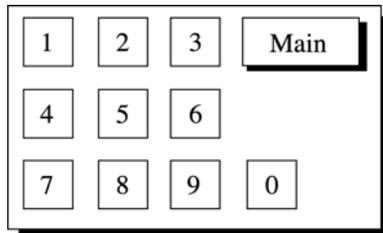

## Clinician Menu

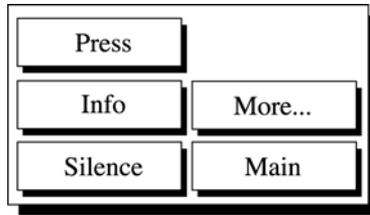

Press. Selection of this button displays a panel for setting the default BP target inflation pressure. Adjusting the default target pressure will automatically update the current inflation target pressure and will be used for the next reading. The range of adjustment is 100 mmHg to 180 mmHg, and the setting is retained when the Monitor is turned off.


The initial target pressure can be set between 100 and 180 mmHg in 5 mmHg increments. The factory default is 150 mmHg for adults and 110 for neonates. This is indicated by the "AUTO" label at the end of the adjustable range. When the target pressure is changed, the next determination will use the new target inflation value if no systolic is available. When adjusted under the Clinician menu of the Service mode, the target pressure is adjusted permanently.

Retain Values. The Retain Values setting indicates the amount of time that the manual BP and Temp measurement values stay on screen. This is also the setting for the amount of time that Spot Mode values are displayed after a freeze period.

OK. Selection of this button returns the user to the Clinician menu.


Info. Selection of this button causes the most recent calibration date of the NIBP system to be displayed. Selection of **OK** returns the user to the Clinician menu.

# Using the Menu System

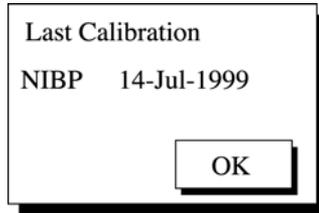

*OK.* Selection of this button returns the user to the Clinician menu.

*More...* Selection of this button displays the More... menu, which allows the user to permanently change default mode settings.

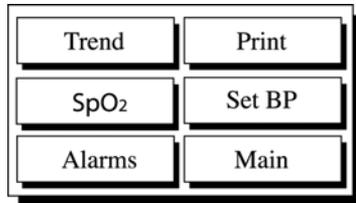

*Trend.* Selection of this button displays the message: **Automatically clear trend on power-up?** Selection of **Yes** overrides the default setting by clearing all trends on power-up and returns the Monitor to the More... menu. Selection of **No** retains the default setting by saving all trends after power-off and returns the Monitor to the More... menu. Selection of **Cancel** returns the user to the More... menu.

*Print.* Selection of this button displays the message: **Restore Print mode on power-up?** Selection of **Yes** restores the Print mode to the default setting (previous user-selected mode) after power-off and returns the Monitor to the More... menu. Selection of **No** restores the Print mode to the Manual mode after power-off and returns the Monitor to the More... menu. Selection of **Cancel** returns the Monitor to the More... menu.

*SpO<sub>2</sub>.* Selection of this button displays the message: **Enter SpO<sub>2</sub> Configuration Mode?** Selection of **No** returns the Monitor to the More... menu. Selection of **Yes** brings up the SpO<sub>2</sub> menu.

For monitors with Nellcor SpO<sub>2</sub> technology, the SpO<sub>2</sub> configuration menu includes Response Mode and *SatSeconds*<sup>TM</sup>. Response Mode can be set to “1” or “2”. The factory default is “2”.

*SatSeconds*<sup>TM</sup> technique limit-controls the time that the % SpO<sub>2</sub> level may fall outside the alarm before an audible alarm sounds. Choose either **0, 10, 25, 50, or 100** seconds. If 0 is chosen this feature is disabled

For monitors with Masimo SpO<sub>2</sub> technology, sensitivity can be set to Normal (the factory default) or Max. FastSAT can be set to On or Off.

*Set BP.* Selection of this button displays the message: **Restore BP mode on power-up?** Selection of **Yes** restores the BP mode to previous user-selected mode after power-off and returns the Monitor to the More... menu. Selection of **No** restores the BP mode to the default setting of Manual after power-off and returns the Monitor to the More... menu. Selection of **Cancel** returns the Monitor to the More... menu.

*Alarms.* Selection of this button displays the message: **Enter alarm configuration mode?** Selection of **No** returns the Monitor to the More... menu. Selection of **Yes** brings up the Alarms menu. Selection of **Reset** changes all alarm limits back to the factory defaults and returns the Monitor to the More... menu. Selection of **Save** permanently saves the user-selected alarm limits and returns the Monitor to the More... menu. Selection of **Cancel** returns the Monitor to the More... menu.


SYS	DIA	MAP	BPM	SpO <sub>2</sub>	Reset
180	130	140	160	OFF	Save
30	15	50	40	90	Cancel

*Main.* Selection of this button returns the user to the Main menu.

# Using the Menu System

Silence. Selection of this button will cause all alarms except the FAILSAFE alarm to be muted. A confirmation menu will appear in Area 2 of the LED. Selection of either **Yes** or **No** returns the user to the Clinician mode menu. If silence is confirmed, the Alarm Silence button (26) illuminates and alarms are permanently muted. If silence is not confirmed, the alarm will be audible.


**Caution: Alarms will be muted until either the Monitor is switched off and on again or the Alarm Silence button (26) is pressed.**

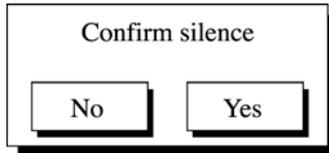

Main. Selection of this button returns the user to the Main menu.

## ***Error and Warning Messages***

The error panel appears in Area 2 of the LCD and indicates the error and its code, if it has one. In this example, a limit violation alarm (which has no error code) has occurred. A list of alarm error messages and their codes is in Appendix B.

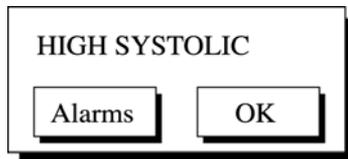

Alarm conditions are addressed in two ways: the Alarms button and OK button.

## **Alarms Button**

Selection of this button takes the user to the Alarms menu, where the alarm limits can be adjusted. This button is available only when a parameter alarm limit has been violated.

### **OK Button**

Selection of this button acknowledges the error. The Monitor clears the identified error and then returns the user to the Main menu.



## **TURBO★TEMP Temperature**

Scale	°Fahrenheit (F) °Celsius (C)
Range	
Predictive mode	Max: 41.1° C; 106.0° F Min: 33.6° C; 96.0° F
Monitor mode	Max: 41.1° C; 106.0° F Min: 26.7° C; 80.0° F
Monitor mode accuracy	± 0.1° C ± 0.2° F (when tested in a calibrated liquid bath; meets ASTM E1112, Table 1, in range specified)
Predictive mode accuracy	± 1.0° F ± 0.6° C
Determination time	approximately 10 seconds

**Use only IVAC probes and probe covers. The size, shape, and thermal characteristics of the probe covers can affect the performance of the instrument. Inaccurate readings or retention problems may occur unless IVAC probes and probe covers are used. Refer to Appendix D for reorder codes.**

### **IVAC® Patents**

U.S. D300,728; D300,909

## **NELLCOR SpO<sub>2</sub>**

### **Measurement Range**

SpO <sub>2</sub>	1 to 100%
Pulse Rate	20 to 250 beats/min
Perfusion Range	0.03 to 20%

### **Accuracy and Motion Tolerance**

#### *Saturation*

Without Motion - Adults*	70 to 100% ±2 digits
Without Motion - Neonate*	70 to 100% ±3 digits
With Motion - Adults/Neo**	70 to 100% ±3 digits
Low Perfusion	70 to 100% ±2 digits

# Appendix A

<i>Pulse Rate</i>	
Without Motion	20 to 250 beats/min $\pm 3$ digits
With Motion	normal physiologic range 55 to 125 beats/min $\pm 5$ digits
Low Perfusion	20 to 250 beats/min $\pm 3$ digits

\*Adult specifications are shown for OXIMAX MAX-A and MAX-N sensors. Neonate specifications are shown for OXIMAX MAX-N. Saturation accuracy will vary by sensor type.

\*\*Applicability: OXIMAX MAX-A, MAX-AL, MAX-P, MAX-I, and MAX-N sensors.

## Default Settings

SpO <sub>2</sub> (%)	HIGH: 100
SpO <sub>2</sub> (%)	LOW: 90
Response Mode	2 (for Mode 2: Fast Response)
Sat Seconds	0
Audible Indicator	Pitch changes continuously with saturation; volume from 0 (off) to 9
Waveforms	Pulse plethysmograph waveform on LCD gain compensated
Sensor Connect/ Disconnect From Patient	Monitor detect attachment or disconnection of sensor from patient within 15 s
Pulse Detection	Monitor will detect pulse or enter no signal state within 15 s of being attached to patient
Loss of Pulse	Monitor will detect loss of pulse from patient and enter no signal state within 10 s
<u>Sensor Light Source</u>	
Wavelength	Infrared: 890 nm (nominal) Red: 660 nm (nominal)
Power Dissipation	52.5 mW (max)

## Nellcor Patents

4,621,643; 4,653,498; 4,700,708; 4,770,179; 4,802,486;  
4,869,254; 4,928,692; 4,934,372; 5,078,136; 5,351,685;  
5,421,329; 5,485,847; 5,533,507; 5,577,500; 5,803,910;  
5,853,364; 5,865,736; 6,083,172; Re. 35,122 and,  
international equivalents.

## Masimo SET SpO<sub>2</sub>

### Measurement Range

SpO <sub>2</sub>	1 to 100%
Pulse Rate	25 to 240 beats/min
Perfusion Range	0.02 to 20%

### Accuracy and Motion Tolerance

#### *Saturation*

Without Motion - Adult/Ped*	70 to 100% ±2 digits
Without Motion - Neonate*	70 to 100% ±3 digits
With Motion - Adult/Ped/Neo**†	70 to 100% ±3 digits
Low Perfusion‡	70 to 100% ±2 digits 0 to 69% unspecified

#### *Pulse Rate*

Without Motion	25 to 240 beats/min ±3 digits
With Motion	normal physiologic range 25 to 240 beats/min ±5 digits

\*The Masimo SET<sup>®</sup> SpO<sub>2</sub> parameter with LNOP-Adt sensors has been validated for no motion accuracy in human blood studies on healthy adult volunteers in induced hypoxia studies in the range of 70-100% SpO<sub>2</sub> against a laboratory co-oximeter and ECG monitor. This variation equals plus or minus one standard deviation. Plus or minus one standard deviation encompasses 68% of the population.

\*\*The Masimo SET<sup>®</sup> SpO<sub>2</sub> parameter with LNOP-Adt sensors has been validated for motion accuracy in human blood studies on healthy adult volunteers in induced hypoxia studies while performing rubbing and tapping motions at 2 to 4 Hz at an amplitude of 1 to 2 cm and a non repetitive motion before 1 to 5 Hz at an amplitude of 2 to 3 cm in induced hypoxia studies in the range of 70-100% SpO<sub>2</sub> against a laboratory co-oximeter and ECG monitor. This variation equals plus or minus one standard deviation. Plus or minus one standard deviation encompasses 68% of the population.

†The Masimo SET<sup>®</sup> SpO<sub>2</sub> parameter with LNOP-Neo Pt sensors has been validated for neonatal motion accuracy in human blood studies on neonates while moving the neonate's foot at 2 to 4 cm against a laboratory co-oximeter and ECG monitor. This validation equals plus or minus, one

# Appendix A

standard deviation. Plus or minus one standard deviation encompasses 68% of the population.

‡The Masimo SET® SpO<sub>2</sub> parameter has been validated for low perfusion accuracy in bench top testing against a Biotek Index 2 simulator and Masimo's simulator with signal strengths of greater than 0.02% and a % transmission of greater than 5% for saturations ranging from 70 to 100%. This variation equals plus or minus, one standard deviation. Plus or minus one standard deviation encompasses 68% of the population.

## *Masimo® Sensor Accuracy*

Sensor Model SpO<sub>2</sub> Range 70% - 100%

### *LNOP*

LNOP-ADT	±2 digits
LNOP-ADT Long	±2 digits
LNOP-PDT	±2 digits
LNOP-NEO	±3 digits
LNOP-NEO PT	±3 digits
LNOP-DCI (reusable)	±2 digits
LNOP-DCSC (reusable)	±2 digits
LNOP-DCIP (reusable)	±2 digits
NRI25 (reusable)	±2 digits

### *Resolution*

Saturation (% SpO <sub>2</sub> )	1%
Pulse Rate (bpm)	1

### *Low Perfusion Performance*

>0.02% Pulse Amplitude and % Transmission >5%	Saturation (% SpO <sub>2</sub> ) ±2 digits Pulse Rate ±3 digits
--	--

## **Interfering Substances**

Carboxyhemoglobin may erroneously increase readings. The level of increase is approximately equal to the amount of carboxyhemoglobin present. Dyes, or any substance containing dyes, that change usual arterial pigmentation may cause erroneous readings.

## **Sensor Light Source**

Wavelength	Infrared: 905 nm (nominal) Red: 660 nm (nominal)
Power Dissipation	Infrared: 22.5 mW (max) Red: 27.5 mW (max)

### **Default Settings**

SpO <sub>2</sub> (%)	HIGH: 100
SpO <sub>2</sub> (%)	LOW: 90
Sensitivity Mode	2 (for low perfusion-Default)
Averaging Time	12 seconds
FastSAT Mode	0 (for Off)

### **Masimo Patents**

5,482,036; 5,490,505; 5,632,272; 5,685,299; 5,758,644;  
5,769,785; 6,002,952; 6,036,642; 6,067,462; 6,206,830;  
6,157,850, and international equivalents.

### **Mechanical**

Dimensions	Height: 9.8 in (25.0 cm) Width: 9.8 in (24.8 cm) Depth: 6.9 in (17.5 cm)
Weight, Including Battery	7.8 lb (3.5 kg)
Mountings	Self-supporting on rubber feet or pole mountable
Portability	Carried by recessed handle or pole mounted
Classification Information	Mode of operation: continuous Degree of protection against harmful ingress of water: Drip- proof IPX1

### **Power Requirements**

**AC input voltage:** 100-240  
VAC, 50 / 60 Hz (nominal)  
90 ~ 253 VAC, 47 ~ 63 Hz  
(range), 50VA.

**Protection against electrical  
shock:** Class 1

**DC input voltage:** 24 VDC  
(nominal), 12-30 VDC, 36VA,

# Appendix A

supplied from a source conforming to IEC 601-1. AC input is protected by two internal fuses, replaceable by qualified service personnel only. DC input line is protected by an internal auto-resetting fuse.

**Battery:** 12 volt, 2.3 amp-hours protected by internal auto-resetting fuse.

**Minimum operation time:** 2 hrs (5 min cycle with adult cuff at 25 °C with power save mode enabled) from full charge.

**Time for full recharge:** 1 hr 50 min from full discharge when the Monitor is switched off and 8 hrs when the Monitor is switched on.

## Environmental

Operating Temperature

+ 5 °C to + 40 °C  
(+ 41 °F to + 104 °F)

Operating Atmospheric Pressure

700 hPa to 1060 hPa

Storage Temperature

- 20 °C to + 50 °C  
(- 4 °F to + 122 °F)

Storage/Transportation Atmospheric Pressure

500 hPa to 1060 hPa

Humidity Range

0% to 95% noncondensing

Radio Frequency

Complies with IEC Publication 601-1-2 (April 1993) Medical Electrical Equipment, Electromagnetic

Compatibility Requirements  
and Tests and CISPR 11  
(Group 1, Class A) for radiated  
and conducted emissions.

## **IPX1**

The DINAMAP<sup>®</sup> PRO Monitor is protected against vertically falling drops of water and conforms to IEC-529 standard at level of IPX1. Vertically falling drops of water shall have no harmful.

# Appendix B

## ***Alarm Codes***

All alarm indications are accompanied by an audible signal unless Alarm Silence is selected.

A microprocessor system failure will generate a high-pitched audible alarm regardless of the setting of the Alarm Silence switch.

There are three categories of alarms: patient alarms, system alarms, and failsafe alarm.

## ***Patient Alarms***

Patient alarms include those alarms issued when the patient's systolic pressure, diastolic pressure, pulse rate, or oxygen saturation is outside the set limits. Whenever one of these conditions occurs, the associated display (SYSTOLIC, MAP, DIASTOLIC, PULSE, or SpO<sub>2</sub>) will flash the most recent reading and an audible alarm will be issued.

Pressing the Alarm Silence switch (causing the integral LED to be lit) silences the audible alarm for 2 minutes, but the alarm display reading and SILENCE LED indicator will continue to flash at the same rate.

## ***System Alarms***

System alarms alert the operator to certain abnormal conditions or internal system failures. Pressing the rotor cancels the alarm information box which is displayed on the LCD. Codes for different procedural and system alarms are on the next page.

## ***Failsafe Alarm***

The failsafe alarm, which is the most powerful alarm of the PRO Monitor, indicates a serious failure of the Monitor. This alarm occurs immediately upon any failure of a self-test and indicates system failure. When the failsafe alarm occurs, the Monitor disables all features to ensure patient safety.

### ***Hierarchy of Alarms***

Alarms in the DINAMAP<sup>®</sup> PRO Monitor are in three priority levels. They are:

<b><u>Alarm</u></b>	<b><u>Priority Level</u></b>
Failsafe	1
Patient and system	2 (High priority alarm)
Low battery	3

The Priority 1 alarm (i.e., Failsafe) will override any other alarm. Priority 2 alarms will override only the low battery alarm. The low battery alarm will not override any other alarm.

# Appendix B

Procedural and Error Alarm Codes							
Alarm Code	LED Display	LCD Description	Audible Tone and Volume	Effect of Alarm Silence Switch	Effect of Clear via SelectKnob	Probable Cause	
N99	No change	N99 - NIBP FAILED	High priority alarm. Volume adjustable	2 minutes silence	Clear	Unable to make an NIBP determination due to insufficient signal	
N55	No change	N55 - TIMEOUT: PRESS	High priority alarm. Volume adjustable	2 minutes silence	Clear	One cuff pressure for > 1 minute. Motion artifact	
N44	No change	N44 - TIMEOUT: TOTAL	High priority alarm. Volume adjustable	2 minutes silence	Clear	Determination time > 2 minutes. Motion artifact	
N33	No change	N33 - TIMEOUT: INFLT	High priority alarm. Volume adjustable	2 minutes silence	Clear	Inflation time > 40 seconds or air leak detected	
N00	No change	N00 - OVER PRESSURE	High priority alarm. Volume adjustable	2 minutes silence	Clear	Overpressure detected	

Procedural and Error Alarm Codes (cont.)							
Alarm Code	LED Display	LCD Description	Audible Tone and Volume	Effect of Alarm Silence Switch	Effect of Clear via SelectKnob	Probable Cause	
P55	No change	P55 - SpO <sub>2</sub> NO SIGNAL	High priority alarm. Volume adjustable	2 minutes silence	Clear	No or very low SpO <sub>2</sub> signal. Check or reposition sensor	
P00	No change	P00 - SpO <sub>2</sub> SENSOR	High priority alarm. Volume adjustable	2 minutes silence	Clear	SpO <sub>2</sub> sensor not connected. No sensor code detected. Sensor failure	
No Code	No change	SpO <sub>2</sub> PLACEMENT? (Masimo Only)	High priority alarm. Volume adjustable	2 minutes silence	Clear	SpO <sub>2</sub> signal weak or noisy. Sensor failure	
No Code	Values zeroed	SpO <sub>2</sub> SENSOR OFF?	High priority alarm. Volume adjustable	2 minutes silence	Clear	SpO <sub>2</sub> Sensor off patient	
E33	No change	E33 - TEMP: FAIL	High priority alarm. Volume adjustable	2 minutes silence	Clear	Temperature probe not connected or inoperable	
E00	No change	E00 - TEMP: FAIL	High priority alarm. Volume adjustable	2 minutes silence	Clear	Predictive temperature determination > 60 sec or attempting axillary temp	
No Code	Blank	HIGH TEMP	High priority alarm. Volume adjustable	2 minutes silence	Clear	Predictive temperature exceeds upper range	
No Code	No change	TEMP PROBE BROKEN	High priority alarm. Volume adjustable	2 minutes silence	Clear	Broken temperature probe	

# Appendix B

Procedural and Error Alarm Codes						
Alarm Code	LED Display	LCD Description	Audible Tone and Volume	Effect of Alarm Silence Switch	Effect of Clear via SelectKnob	Probable Cause
No Code	No change	LOW BATTERY, Flashing battery icon	3 beeps every 10 seconds, adjustable volume	2 minutes silence	No effect	Replace or recharge battery. From onset of alarm. 5 NIBP measurements available. Beep rate increases linearly as battery discharges
No Code	Blank	LOW BATTERY - SYSTEM DISABLED	Steady tone, maximum volume	No effect	No effect	Replace or recharge battery. NIBP measurement disabled
No Code	No change	PRINTER - NO PAPER	High priority alarm. Volume adjustable	2 minutes silence	Clear	Paper ran out or printer door open
No Code	Values posted	NIBP RANGE ERROR	High priority alarm. Volume adjustable	2 minutes silence	Clear	NIBP algorithm returned value outside specified accuracy range
Other: N, P, E, I, S	Blank	Error code, description	Steady tone, maximum volume	No effect	No effect	Internal system fault




# Appendix C

## ***Principles of Noninvasive Blood Pressure Determination***

The oscillometric method of determining NIBP is accomplished by a sensitive transducer which measures cuff pressure and minute pressure oscillations within the cuff. A single determination (in normal mode) is initiated before taking repeated determinations in auto or stat mode. As a determination is taken, the algorithm stores the pattern of the patient's oscillation size as a function of the pressure steps. In auto and stat mode, as few as six pressure steps may be necessary to complete the determination process. When employing fewer pressure steps, the system uses the stored information from the previous blood pressure determination to decide the best pressure steps to take. The algorithm measures the consistency of pulse size to tell if the oscillations taken at a step are good and if more steps are needed.

The first determination initially pumps up to a target pressure of about 150 mmHg for adult patients depending on initial target pressure preset. After inflating the cuff, the NIBP parameter begins to deflate, the oscillations versus cuff pressure is measured, and, finally, systolic pressure, mean pressure, and diastolic pressure are calculated and the screen is updated.

During an NIBP determination, the parameter deflates the cuff one step each time it detects two pulsations of relatively equal amplitude. The time between deflation steps depends on the frequency of these matched pulses (pulse rate of the patient). However, if the monitor is unable to find any pulse within several seconds, it will deflate to the next step. The process of finding two matched pulses at each step provides artifact rejection due to patient movement and greatly enhances the accuracy of the monitor. The figure shows a full determination sequence for an adult patient. In stat mode, the artifact rejection technique of two matching pulses at each step is disabled and some steps may have only one pulse.

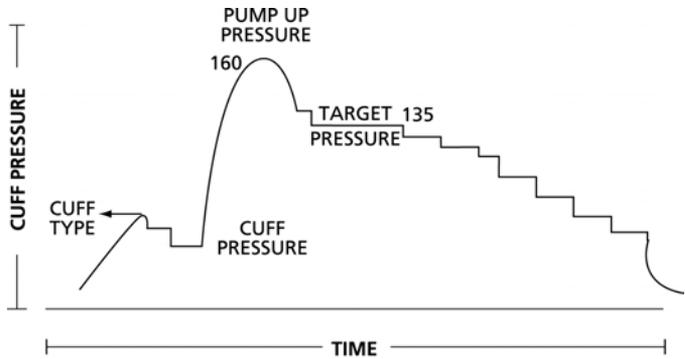


### Full NIBP Determination Sequence for Adult

At each step the microprocessor stores cuff pressure, the matched pulse amplitude, and the time between successive pulses. The stepped deflation and matched pulse detection continues until diastolic pressure is determined or total cuff pressure falls below 8 mmHg. The parameter then deflates the cuff (to zero detected pressure), analyzes the stored data, and updates the screen.

The operating cycle is composed of four parts: inflation time, deflation time, evaluation time, and wait time. Wait time, which varies from mode to mode, is affected by the cycle time (auto mode) or operator intervention (manual mode). The figure shows the basic operating cycle for an NIBP determination.

# Appendix C

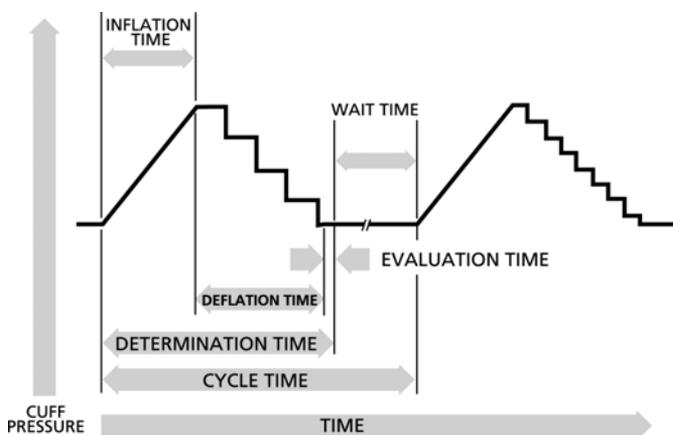

## NIBP - Auto Mode

### Systolic Search

If systolic pressure is not found, the NIBP parameter can search at cuff pressures higher than the initial target pressure. The parameter will inflate the cuff above the initial target pressure to get better data in the systolic region. The maximum pressure allowed in systolic search is limited by the normal range for cuff pressures. In any operating mode, if a patient's systolic pressure exceeds the inflation pressure, the parameter will begin normal deflation sequence, detect the absence of a systolic value, stop deflation, reinflate to a higher (than initial) inflation pressure (290 mmHg maximum), and resume the normal deflation sequence.

Do not use the auscultatory method to verify the accuracy of the NIBP parameter. Auscultatory method (using cuff and stethoscope) calculates the mean pressure value from audible sounds at systolic and diastolic, but the NIBP parameter method detects all three values.

### Reverting and Accelerated Determination

When the determination data does not agree with the previous determination, there may be a reversion to a determination with more steps. The algorithm will try to

make an accelerated determination of blood pressure if it has been 16 minutes or less since the last determination.

# Appendix D

## Reorder Codes

### **PRODUCT**

	<b>PROD CODE</b>	<b>FCAT</b>
DINAMAP <i>PRO</i> Monitor Operations Manual-English	2009802-001	
DINAMAP <i>PRO</i> Monitor Service Manual	2009829-001	
Battery, 12V Lead Acid	2010422-001	
Printer Paper (box of 10)	089100	E9050KP
DINAMAP Rolling Stand	003215	E9050JB
Nurse Call	487208CR	

### **NIBP:**

Air Hose 12 ft Adult/Pediatric, Screw Connector	107365	E9050LH
Air Hose 24 ft Adult/Pediatric, Screw Connector	107366	E9050LJ
Air Hose 12 ft Neonatal, Quick Disconnect	107368	E9050LK
Air Hose 12 ft Adult/Pediatric, Quick Disconnect	88847	E9050KN

### **CUFF Assortment Packs:**

CLASSIC-CUF <sup>®</sup> Assortment Pack Includes: 1 each: Infant, Child, Small Adult, Adult, Large Adult, Thigh Cuff	2692	E2692J
---	------	--------

CLASSIC-CUF <sup>®</sup> Assortment Pack, Neonate Includes: 2 Neo #1, 3 Neo #2, 5 Neo #3, 5 Neo #4, 5 Neo #5	2693	E2693J
--	------	--------

SOFT-CUF <sup>®</sup> Assortment Packs: Includes: 1 Infant, 1 Child, 2 Small Adult, 2 Adult, 2 Large Adult, 1 Thigh, 1 Adult Long Cuff	2695	E2695J
--	------	--------

SOFT-CUF <sup>®</sup> Assortment, Neonate Includes: 2 Neo #1, 3 Neo #2, 5 Neo #3, 5 Neo #4, 5 Neo #5	2694	E2694J
--	------	--------

DURA-CUF <sup>®</sup> Assortment Packs: Includes: 1 each: Infant, Child, Small Adult, Adult, Large Adult, Thigh Cuff	2699	E2699J
--	------	--------

DURA-CUF <sup>®</sup> Assortment Pack, Adult Includes: 1 each: Infant, Child, Small Adult, Adult, Large Adult, Thigh Cuff	2698	E2698J
---	------	--------

DURA-CUF <sup>®</sup> Assortment Pack, Child Includes: 2 Infant, 3 Child, and 1 Small Adult Cuff	2697	E2697J
---	------	--------

Additional Blood Pressure Cuff Codes are available through: [www.gemedical.com](http://www.gemedical.com)

### **TEMPERATURE:**

IVAC <sup>®</sup> TURBO★TEMP Oral Temp Probe, Long Cord	2008774-001	
IVAC <sup>®</sup> TURBO★TEMP Rectal Temp Probe, Long Cord	2008775-001	
IVAC <sup>®</sup> Temperature Probe Covers	88015	E9050KK

### **SPO2:**

#### **NELLCOR<sup>®</sup>:**

Pulse Oximeter Cable DOC-10	2008773-001	
DuraSensor Adult Oxygen Sensor	DS100A	

#### **Masimo<sup>®</sup>:**

Adult Reusable Sensor, 1/Bx (NR125)	2009745-001	
Cable (PC08)	2009743-001	



# Appendix E

## ***Warranty, Service, and Spare Parts***

**Warning: There are no user serviceable parts inside the DINAMAP® PRO Monitor. Refer all servicing to qualified personnel.**

### ***Warranty***

All repairs on products under warranty must be performed or approved by GE Medical Systems *Information Technologies*. Unauthorized repairs will void the warranty.

Products not covered by warranty should be repaired only by qualified electronics service personnel.

### ***Assistance and Parts***

If the product malfunctions or if assistance, service, or spare parts are required, contact GE Medical Systems *Information Technologies* Technical Support.

Before contacting GE Medical Systems *Information Technologies* it is helpful if you can duplicate the problem and check and confirm the operation of all accessories to ensure that they are not the cause of the problem.

When calling, please have the following information at hand:

- product name and model number and complete description of the problem
- the serial number of your Monitor
- your name and address
- a purchase order number if out-of-warranty repairs or spare parts are required
- your GE Medical Systems *Information Technologies* account number, if applicable
- the 6-digit part number for spare or replacement parts

### ***Repairs***

If your product requires warranty, extended warranty, or non-warranty repair service, call GE Medical Systems *Information Technologies* and a representative will assist you.

Estimates for non-warranty repairs are provided at no charge; however, the product must be sent to GE Medical Systems *Information Technologies* for an estimate.

To facilitate prompt service in cases where the product has external chassis or case damage, please advise the representative when you call.

The representative will record all necessary information and will provide a Return Authorization Number. Prior to returning any product for repair, a Return Authorization Number must be obtained.

### ***Packing Instructions***

If you have to return goods for service, follow these recommended packing instructions:

- Remove all hoses, cables, sensors, power cords, and ancillary products from the Monitor before packing.
- Wherever possible use the original shipping carton and packing materials.
- Observe the environmental conditions detailed in Appendix A.

It is recommended that all returned goods be insured. Claims for loss or damage to the product must be initiated by the sender.

### ***Service Manuals***

Service Manuals containing calibration and repair information can be ordered from GE Medical Systems *Information Technologies*. These manuals also include full schematic diagrams, assembly drawings, and spare parts lists. Refer to Appendix D for the reorder number of the Service Manual.

# Appendix F

## **Maintenance**

### **Cleaning the Monitor**

The Monitor and accessories are to be kept clean and used according to the instructions provided here and in the Service Manual.

The exterior of the Monitor may be wiped clean with a soft cloth slightly dampened with mild detergents. The Monitor and accessories should be inspected once yearly for wear and damage.

- Do **not** immerse unit.
- Do **not** clean with isopropyl alcohol or other solvents.
- Do **not** immerse hoses.

### **Cuff Cleaning and Disinfection**

#### *General*

The cuff must be thoroughly cleaned with the specified detergent before reuse. The additional use of household bleach as described below provides at least intermediate-level disinfection.

- Apply cuff hose plugs before cleaning.
- The following cleansing procedure was repeated 20 times on DURA-CUF® Blood Pressure Cuffs and once on SOFT-CUF® Blood Pressure Cuffs without affecting the performance of the cuff.
- While this procedure is adequate for cleaning/disinfection, it may not remove all stains.
- Do **not** immerse hoses.
- Do **not** immerse cuffs without prior application of cuff hose caps.

#### *Materials*

- Enzymatic detergent such as ENZOL\* enzymatic detergent (US) or Cidezyme\* enzymatic detergent (UK)
- Distilled water
- 10% solution of household bleach (5.25% sodium hypochlorite) in distilled water
- Soft cloths and soft-bristled brushes

\*Trademark

- Spray bottles

#### *Procedure*

1. Prepare the enzymatic detergent according to the manufacturer's instructions and the 10% bleach solution, in separate spray bottles.
2. Spray the detergent liberally on device. If the material is dried on, allow the cuff to sit for 1 minute. For soil on the soft part of the closure or the cuff itself, wipe the material off with a soft cloth. For persistent contamination on the soft part of the closure, use a soft-bristled brush to loosen particles. Rinse with copious amounts of distilled water. Repeat until no visible contamination remains. For soil on the hook part of the closure, use a soft-bristled brush to remove the material, and rinse with copious amounts of distilled water. Repeat until no visible contamination remains.
3. Spray the 10% bleach solution on the affected area until the area is saturated. Allow the cuff to sit for 5 minutes.
4. Wipe away any excess solution and rinse the cuff again with distilled water. Allow 2 hours for drying.

The user has the responsibility to validate any deviations from the recommended method of cleaning and disinfection.

For additional information on infection control procedures, contact GE Medical Systems *Information Technologies* Technical Support.

#### **Temperature Devices**

Do not immerse predictive temperature probes. The probe may be cleaned with an alcohol solution. Use a cloth or sponge—just damp, not wet—and avoid getting any liquid into the interior of the probe.

#### **SpO<sub>2</sub> Sensors**

Adhesive sensors are sterile and for single use only. Reusable sensors should be cleaned before reuse with a 70% alcohol solution. Do not immerse the sensor completely in water, solvents, or cleaning solutions (because the connector is not waterproof). Do not sterilize the sensor by irradiation, steam, or ethylene oxide. If disposable sensors or their packaging

# Appendix F

are damaged, they must be disposed of as advised in Appendix F.

## ***Storage and Battery Care***

If it becomes necessary to store the Monitor for an extended period of time, first fully charge then remove the battery. Then store the Monitor and the battery in the original packaging materials.

Batteries should always be fully charged before being placed in storage. Even after 6 months of storage, a fully charged battery can retain about 80% of its charge. A fully charged battery in good condition will provide sufficient power to operate a Monitor for approximately 2 hours, including temperature and BP measurements made at 5-minute intervals.

It is best to keep the battery charged as fully as practical and never store the Monitor with the battery in a discharged condition. When the battery will no longer hold a charge, remove and replace it with one of the same part number. Failure to replace the battery with the same GE Medical Systems *Information Technologies* part number may result in shorter battery life.

To charge the battery, insert the plug from either the AC mains power cord or the AC-DC power converter into an appropriate AC outlet. The battery will charge regardless of the position of any switches.

Battery charging will take place as long as the Monitor remains connected to an external AC power source. A battery that is fully discharged can be fully recharged in 1 hour 50 minutes when the Monitor is switched off or 8 hours if the Monitor is switched on.

## **Cautions**

- **To ensure that the battery will be ready for portable operation, keep the Monitor connected to a mains supply whenever possible.**

- **Repeated failure to fully charge the battery will result in a significant reduction in battery life.**
- **The expected lifetime of the battery largely depends on the way in which the Monitor is used. If the battery is allowed to completely discharge before being fully recharged, the battery should survive around 200 recharge cycles. If the battery is used in such a way that it never becomes more than one third discharged and is fully recharged whenever possible, it can survive up to 1200 cycles. This means that by thoughtful usage, the lifetime of the battery can be extended up to six times.**

Replacement batteries may be obtained from GE Medical Systems *Information Technologies*.

**Note:** The replacement part number of the battery is 633132. Do not use other types.

### ***Fuses***

The Monitor contains five fuses. Two AC line input fuses are mounted internally and are replaceable only by qualified service personnel. The remaining three fuses are auto-resettable and mounted within the Monitor. These fuses protect the low voltage DC input, the battery, and the +5 V output on the host port connector.

### ***Calibration***

Calibration of the BP parameter should be checked at least once a year or when there is doubt about the validity of the readings.

### ***Leak Testing***

A leak test of the BP parameter should be performed at least once a year or when there is doubt about the validity of the pressure readings.

**Caution: Refer calibration and leak testing to qualified service personnel. Full calibration details are available in the DINAMAP PRO Monitor Service Manual, available from GE Medical Systems *Information Technologies*.**

# Appendix F

## ***Disposal of Product Waste***

As you use the PRO Monitor, you will accumulate solid wastes that require proper disposal or recycling. These include batteries, patient applied parts, and packaging material.

### **Batteries**

#### **Caution: Do not incinerate batteries.**

The sealed, rechargeable backup battery contains lead and can be recycled. The rechargeable memory battery is of the Nickel Metal Hydride form. Discharge this battery prior to disposal. Place the battery in packaging which electrically isolates its contents. Do not puncture or place the battery in a trash compactor. Do not incinerate the battery or expose it to fire or high temperatures. Dispose in accordance with regional body controlled guideline.

### **Patient Applied Parts**

Certain patient applied parts, such as those with adhesive (disposable SpO<sub>2</sub> sensors), are intended for single use and should be disposed of properly as medical waste in accordance with regional body controlled guideline.

Other patient applied parts, such as blood pressure cuffs, should be cleaned according to instructions. Inspect reusable applied parts for wear, replace as necessary, and dispose of used product as medical waste in accordance with regional body controlled guideline.

### **Packaging Material**


Retain original packaging materials for future use in storing or shipping the Monitor and accessories. This recommendation includes corrugated shippers and inserts.

Whenever possible recycle the packaging of accessories and patient applied parts.



# Appendix G

## Host Port Connector (rear panel)

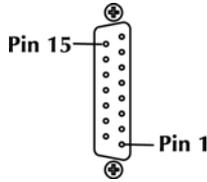

WARNING! Auxiliary equipment connected to the DINAMAP® PRO Monitor will result in the formation of an electromedical system and thus, must comply with the requirements of EN 60601-1-1/ IEC 601-1. All host port signals are NON-ISOLATED and should be connected to equipment conforming to IEC-601-1, configured to comply with IEC 601-1-1 ONLY. Where isolation of data communication is required, GE Medical Systems *Information Technologies* part number ILC232 should be used. If external alarm control is required, GE Medical Systems *Information Technologies* part number 487208 (Isolated Remote Alarm Cable Assembly) should ALWAYS be used. Please refer to the Information Sheet included with the isolated remote alarm cable for details.

**Note:** When using remote alarm, the PRO Monitor should be considered the primary alarm source. The secondary alarm is used for secondary purposes only.

### Pin NumberFunction

1	Ground
2	Inverted TTL Transmit Data
3	Inverted TTL Receive Data
4	Fused +5 volts
5	No connection
6	No connection
7	Ground
8	Remote Alarm
9	RS232 Request to Send (RTS)
10	RS232 Clear to Send (CTS)
11	RS232 Transmit Data (TxD)
12	No connection
13	RS232 Receive Data (RxD)
14	No connection
15	No connection





**European Representative**

GE Medical Systems  
Information Technologies  
GmbH

Postfach 60 02 65

D-79032 Freiburg,

Germany

Tel. +49 761 45 43 - 0

Fax +49 761 45 43 - 233



2009802-001 B